


Manual

Vibratory Sieve Shaker AS 450 control



Translation



Copyright

© Copyright by
Retsch GmbH
Retsch-Allee 1-5
42781 Haan
Germany

Table of Contents

1	Notes on the manual	6
1.1	Explanation of signs and symbols	6
1.2	Disclaimer	6
1.3	Copyright.....	6
2	Safety	7
2.1	Explanations of the Safety Instructions	8
2.2	General Safety Instructions	9
2.3	Repairs.....	10
2.4	Confirmation Form for the Managing Operator.....	11
3	Technical Data	13
3.1	Degree of Protection.....	13
3.2	Emissions.....	13
3.3	Electromagnetic Compatibility (EMC)	14
3.4	Rated Power	14
3.5	Dimensions and Weight.....	14
3.6	Required Floor Space	14
3.7	Receptacle Volume.....	14
3.8	Feed Grain Size	15
3.9	Payload	15
3.10	Suitable Sieve Diameters	15
4	Packaging, Transport and Installation	16
4.1	Packaging	16
4.2	Transport.....	16
4.3	Temperature Fluctuations and Condensation	16
4.4	Conditions for the Installation Site	17
4.5	Electrical Connection	18
4.6	Type Plate Description.....	18
4.7	Removing the Transportation Aid	19
5	First Commissioning	22
5.1	Sieve Clamping Unit "standard"	23
5.2	Sieve Clamping Unit "comfort"	23
5.3	Connecting the Control Unit.....	25
5.4	Wall Mounting of the Control Unit.....	26
6	Operating the Device	28
6.1	Use of the Device for the Intended Purpose.....	28
6.2	Principle of Operation	29
6.3	Views of the device	30
6.3.1	Front.....	30
6.3.2	Back	31
6.4	Switching On / Off.....	31
6.5	Selection of the Test Sieves	32
6.6	Performing a Sieving	32
7	Controlling the Device	33
7.1	Operating Controls, Displays and Functions	33
7.1.1	Start Process	33
7.1.2	Stop Process.....	34
7.2	Pause Process.....	34
7.3	Amplitude	34
7.3.1	Amplitudes in Dependence on the Load.....	35
7.4	Time	36
7.5	Optimisation of Time and Amplitude.....	36
7.6	Interval	37

7.6.1	Interval Time	37
7.7	Programme Mode	38
7.7.1	Select a Program	38
7.7.2	Edit a Program	38
7.7.3	Save a Programme	38
7.8	Signal Tone	39
7.8.1	Operating Hours	39
7.8.2	Software Version	39
8	Wet Sieving	41
8.1	Performing the Wet Sieving	41
9	EasySieve®	44
10	Error Messages and Information Notes	45
10.1	Error Messages	45
10.2	Information Notes	45
11	Return for Service and Maintenance	46
12	Cleaning, Wear and Maintenance	47
12.1	Cleaning	47
12.1.1	Cleaning of Test Sieves	47
12.1.1.1	Cleaning of Test Sieves with Mesh Sizes > 500 µm	47
12.1.1.2	Cleaning of Test Sieves with Mesh Sizes < 500 µm	47
12.1.1.3	Drying of Test Sieves	48
12.2	Wear	48
13	Maintenance	48
13.1.1	Replacing the Fuses	49
14	Accessories	50
14.1	Test Sieves	50
14.1.1	Certificate	50
14.1.2	Calibration Service	51
14.2	Sieving Aids	51
15	Disposal	52
16	Index	53



1 Notes on the manual

This manual provides technical guidelines for the safe operation of the device. Read this manual through carefully before installing, putting into service and operating the device. Reading and understanding this manual is essential for handling the device safely and as intended.

This manual does not contain any repair instructions. Please contact your supplier or contact Retsch GmbH directly if anything is unclear or you have questions about these guidelines or the device, or in the case of any faults or necessary repairs.

You can find further information about your device at <https://www.retsch.com> on the pages for the specific device concerned.

Amendment status:

The document amendment 0011 of the "Vibratory Sieve Shaker AS 450 control" manual has been prepared in accordance with the Machinery Directive 2006/42/EC.

1.1 Explanation of signs and symbols

In this document the following **signs and symbols** are being used:

ⓘ	Reference to a recommendation and/or an important information
→	Reference to a chapter, table or figure
⇒	Action instruction
Name	Software menu function
[Name]	Software button
<Name>	Software checkbox

1.2 Disclaimer

This manual has been prepared with great care. We reserve the right to make technical changes. We assume no liability for personal injuries resulting from the failure to follow the safety information and warnings in this manual. No liability will be assumed for damage to property resulting from the failure to follow the information in this manual.

1.3 Copyright

This document or parts of it or its content may not be reproduced, distributed, edited or copied in any form without prior written permission of Retsch GmbH. Damage claims shall be asserted in the case of infringements.

2 Safety

Safety Officer

The operating company itself must ensure the following with respect to persons authorised to work on the device:

- that they have read and understood all regulations contained in the chapter on safety;
- that they are aware before they start work of all instructions and regulations for the target group related to the work;
- that they have easy access to the manual for this device at all times;
- that they have been familiarised with the safe and correct handling of the device before starting work on it, by means of a verbal introduction by a competent person and/or using this manual.

⚠ Improper operation can lead to personal injuries. The operating company itself is responsible for its safety and that of its staff. The operating company itself must ensure that no unauthorised persons have access to the device.

Target group

All those operating, cleaning or working with or on the device.

This device is a modern, powerful product from Retsch GmbH and has been developed in line with the state-of-the art. The device is safe to use when operated correctly and when following the instructions in this manual.

⚠ People under the influence of intoxicating substances (medications, drugs, alcohol) or who are overtired may not operate the device or work on the device.

2.1 Explanations of the Safety Instructions

The following **warnings** in this manual warn of possible risks and damage:

⚠ DANGER	<small>D1.0000</small>
<p>Risk of fatal injuries Source of danger</p> <ul style="list-style-type: none"> – Possible consequences if the danger is ignored. • Instructions and information on how to avoid the risk. 	

Fatal or serious injuries may result if the “Danger” sign is disregarded. There is a **very high risk** of a life-threatening accident or lasting personal injury. The signal word **⚠ DANGER** is additionally used in the running text or in instructions.

⚠ WARNING	<small>W1.0000</small>
<p>Risk of life-threatening or serious injuries Source of danger</p> <ul style="list-style-type: none"> – Possible consequences if the danger is ignored. • Instructions and information on how to avoid the risk. 	

Life-threatening or serious injuries may result if the “Warning” sign is disregarded. There is an **increased risk** of a serious accident or of a possibly fatal personal injury. The signal word **⚠ WARNING** is additionally used in the running text or in instructions.

⚠ CAUTION	<small>C1.0000</small>
<p>Risk of injuries Source of danger</p> <ul style="list-style-type: none"> – Possible consequences if the danger is ignored. • Instructions and information on how to avoid the risk. 	

Average to slight injuries may result if the “Caution” sign is disregarded. There is an average or slight risk of an accident or personal injury. The signal word **⚠ CAUTION** is additionally used in the running text or in instructions.

NOTICE

N1.0000

Type of damage to property

Source of the damage to property

- Possible consequences if the information is ignored.
- **Instructions and information on how to avoid the damage to property.**

Damage to property may result if the information is disregarded. The signal word **NOTICE** is additionally used in the running text or in instructions.

2.2 General Safety Instructions

CAUTION

C2.0002

Risk of injury

Lack of knowledge of the manual

- The manual contains all safety-related information. Disregarding the manual can therefore lead to injuries.
- **Read the manual carefully before operating the device.**



CAUTION

C3.0015

Risk of injury

Improper modifications to the device

- Improper modifications to the device can result in injuries.
- **Do not make any unauthorised changes to the device.**
- **Only use the spare parts and accessories approved by Retsch GmbH!**

NOTICE

N2.0012

Changes to the device

Improper modifications

- The conformity declared by Retsch GmbH with the European Directives will lose its validity.
- Any warranty claims will be terminated.
- **Do not make any modification to the device.**
- **Use spare parts and accessories that have been approved by Retsch GmbH exclusively.**



2.3 Repairs

This manual does not contain any repair instructions. For safety reasons, repairs may only be carried out by Retsch GmbH or an authorised representative or by qualified service technicians.

In case of repair, please inform...

- ...the Retsch GmbH representative in your country,
- ...your supplier, or
- ...Retsch GmbH directly.

Service address:

2.4 Confirmation Form for the Managing Operator

This manual contains essential instructions for operating and maintaining the device which must be strictly observed. It is essential that they be read by the user and by the qualified staff responsible for the device before the device is commissioned. This manual must be available and accessible at the place of use at all times.

The user of the device herewith confirms to the managing operator (owner) that he has received sufficient instructions about the operation and maintenance of the system. The user has received the manual, has read and taken note of its contents and consequently has all the information required for safe operation and is sufficiently familiar with the device.

The managing operator should for legal protection have the user confirm the instruction about the operation of the device.

I have read and taken note of the contents of all chapters in this manual as well as all safety instructions and warnings.

User

Surname, first name (block letters)

Position in the company

Place, date and signature

Managing operator or service technician

Surname, first name (block letters)

Position in the company

Place, date and signature

3 Technical Data

3.1 Degree of Protection

- IP54 (device)
- IP20 (control unit)

3.2 Emissions

⚠ CAUTION

C4.0011

Possibility of acoustic signals not being heard

Loud sieving noises

- Possible acoustic alarms and voice communication might not be heard.
- **Consider the volume of the sieving noise in relation to other acoustic signals in the work environment. Additional visual signals may be used.**

⚠ CAUTION

C5.0017

Hearing damage

A high sound level may be generated depending on the type of material, the number of sieves, the sieving aid used, the amplitude set and the duration of the sieving



- Excessive noise in terms of level and duration can cause impairments or permanent damage to hearing.
- **Ensure suitable noise protection measures are taken or wear ear protection.**

Sound parameters:

The sound parameters are also influenced by the set amplitude, the number of test sieves and the properties of the sample material.

Example 1:

Number of test sieves:	5
Amplitude:	1.1 mm
Feed material:	Quartz sand (< 1 mm)

At these operating conditions, the workplace related equivalent continuous sound level $L_{eq} = 69 \text{ dB(A)}$.

Example 2:

Number of test sieves:	5
Amplitude:	2.2 mm
Feed material:	Quartz sand (< 1 mm)

At these operating conditions, the workplace related equivalent continuous sound level $L_{eq} = 79 \text{ dB(A)}$.

3.3 Electromagnetic Compatibility (EMC)

- EMC class according to DIN EN 55011: B

Strong electromagnetic interference fields, such as high-power radio transmitters, can have an adverse influence on the amplitude control of the AS 450 control. Once the source of the interference is eliminated, the AS 450 control will return to normal operation by itself.

3.4 Rated Power

~ 1 500 W (VA)

3.5 Dimensions and Weight

- Height without sieve clamping unit: 440 mm
- Height with sieve clamping unit: 1 450 mm
- Width: 705 mm
- Width with clamping unit: 900 mm
- Depth: 635 mm
- Weight without sieve stack, without clamping unit: ~ 220 kg

3.6 Required Floor Space

- Width of the base: 800 mm
- Depth of the base: 700 mm
- No safety clearances required

Location requirements:

The device must be placed on a vibration-free, plane, stable and free surface to avoid transmission of vibrations. A level base ensures the uniform distribution of the sample over the sieve mesh fabric, as well as the stability of the device.

3.7 Receptacle Volume

The maximum receptacle volume (the maximum feed quantity) depends on various factors such as number and aperture size of the test sieves, maximum grain size and width of distribution of the sample material.

Examples for the maximum feed quantity according to DIN 66165 for test sieves of 450 mm (18") in diameter are listed in the following table:

Mesh size	Max. feed quantity	Max. permitted oversize material according to DIN 66165
25 µm	64 cm ³	32 cm ³
45 µm	95 cm ³	48 cm ³
63 µm	127 cm ³	64 cm ³
125 µm	191 cm ³	95 cm ³
250 µm	286 cm ³	143 cm ³
500 µm	445 cm ³	223 cm ³
1 mm	636 cm ³	318 cm ³
2 mm	1 113 cm ³	557 cm ³
4 mm	1 749 cm ³	875 cm ³
8 mm	2 863 cm ³	1 431 cm ³

3.8 Feed Grain Size

Traditional dry sieving is performed in the particle size range of 40 µm to 125 mm. By means of sieving aids or with wet sieving the measurement range can be extended to 20 µm. The maximum feed grain size depends on the sample material, the number and aperture size of the test sieves and the type of the sieving machine.

Examples for the maximum feed grain size according to DIN 66165 are listed in the following table:

Mesh size	Max. feed grain size according to DIN 66165	Mesh size	Max. feed grain size according to DIN 66165
22 µm	710 µm	4 mm	25 mm
45 µm	1 mm	8 mm	45 mm
63 µm	1.4 mm	16 mm	71 mm
125 µm	2.5 mm	22.4 mm	90 mm
250 µm	4 mm	45 mm	150 mm
500 µm	6 mm	63 mm	180 mm
1 mm	10 mm	90 mm	230 mm
2 mm	16 mm	125 mm	300 mm

The Vibratory Sieve Shaker AS 450 control is designed for the measurement range of 25 µm to 125 mm.

3.9 Payload

- Maximum sample quantity: 25 kg
- Maximum sieve stack weight: 50 kg
- Maximum payload: 75 kg (sample material plus test sieves)
- Maximum sieve stack height: 963 mm
- Maximum number of fractions: 9 (test sieves and collecting pan: 450 mm) / 13 (test sieves and collecting pan: 400 mm)

3.10 Suitable Sieve Diameters

- Suitable sieve diameters: 400 mm / 450 mm (18")

4 Packaging, Transport and Installation

4.1 Packaging

The packaging has been adapted to the mode of transport. It complies with the generally applicable packaging guidelines.

NOTICE

N3.0001

Complaint or return

Keeping the packaging

- Inadequate packaging and insufficient securing of the device can jeopardise the warranty claim in the event of a complaint or return.
- **Keep the packaging for the duration of the warranty period.**

4.2 Transport

NOTICE

N4.0017

Damage to components

Transport

- Mechanical or electronic components may be damaged during transport. The device must not be knocked, shaken or thrown during transport.
- **Move the device gently during transport.**

NOTICE

N5.0014

Complaints

Incomplete delivery or transport damage

- The forwarding agent and Retsch GmbH must be notified immediately in the event of transport damage. It is otherwise possible that subsequent complaints will not be recognised.
- **Please check the delivery on receipt of the device for its completeness and intactness.**
- **Notify your forwarding agent and Retsch GmbH within 24 hours.**

4.3 Temperature Fluctuations and Condensation

NOTICE

N6.0016

Damaged components due to condensation

Temperature fluctuations

- The device may be exposed to substantial fluctuations in temperature during transport. The ensuing condensation can damage electronic components.
- **Wait until the device has acclimatised before putting it into service.**

Temporary storage:

Also in case of an interim storage the device must be stored dry and within the specified ambient temperature range.

4.4 Conditions for the Installation Site

NOTICE

N7.0021

Ambient temperature

Temperatures outside the permitted range

- Electronic and mechanical components may be damaged.
- The performance data alter to an unknown extent.
- **Do not exceed or fall below the permitted temperature range (5 °C to 40 °C ambient temperature) of the device.**
- Installation height: max. 2 000 m above sea level
- Ambient temperature: 5 °C – 40 °C
- Maximum relative humidity < 80 % (at ambient temperatures ≤ 31 °C)

For ambient temperatures U_T between 31 °C and 40 °C, the maximum relative humidity value L_F linearly decreases according to $L_F = -(U_T - 55) / 0.3$:

Ambient temperature	Max. rel. humidity
≤ 31 °C	80 %
33 °C	73.3 %
35 °C	66.7 %
37 °C	60 %
39 °C	53.3 %
40 °C	50 %

NOTICE

N8.0015

Humidity

High relative humidity

- Electronic and mechanical components may be damaged.
- The performance data alter to an unknown extent.
- **The relative humidity in the vicinity of the device should be kept as low as possible.**


4.5 Electrical Connection

⚠
WARNING

W2.0015

Risk to life caused by an electric shock
 Connection to socket without a protective earth conductor

- Connecting the device to sockets without a protective earth conductor can lead to life-threatening injuries caused by an electric shock.
- **Always operate the device using sockets with a protective earth conductor (PE).**



N
NOTICE

N9.0022

Electrical connection
 Failure to observe the values on the type plate

- Electronic and mechanical components may be damaged.
- **Connect the device only to a mains supply matching the values on the type plate.**

- ⚠ WARNING** When connecting the power cable to the mains supply, use an external fuse that complies with the regulations applicable to the place of installation.
- Check the type plate for details on the necessary voltage, frequency, and maximum external current source fuse for the device.
 - The listed values must agree with the existing mains supply.
 - Only use the supplied power cable to connect the device to the mains supply.

4.6 Type Plate Description

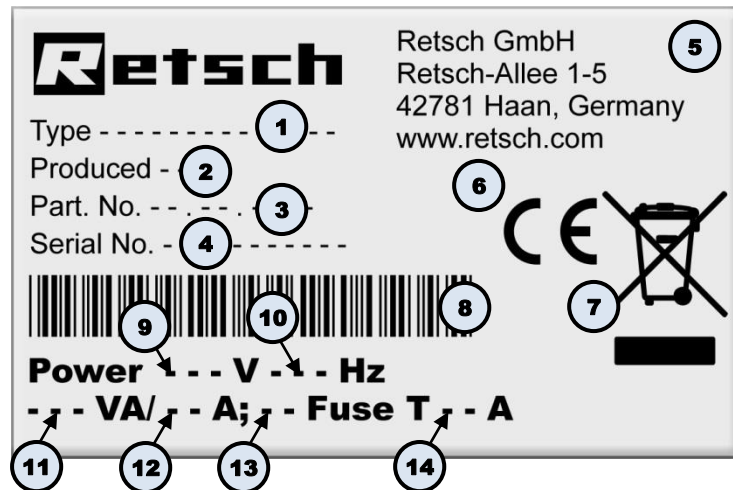


Fig. 1: Type plate

- 1 Device designation
- 2 Year of production
- 3 Part number
- 4 Serial number
- 5 Manufacturer's address
- 6 CE marking
- 7 Disposal label
- 8 Bar code

- 9 Power version
- 10 Mains frequency
- 11 Capacity
- 12 Amperage
- 13 Number of fuses
- 14 Fuse type and fuse strength

⑩ In the case of queries please provide the device designation (1) or part number (3), as well as the serial number (4) of the device.

4.7 Removing the Transportation Aid

⚠ WARNING W3.0005

Risk of injury due to the device falling down
Lifting the device above head height

- The device can fall causing serious injuries when lifted above head height.
- **Never lift the device above head height!**

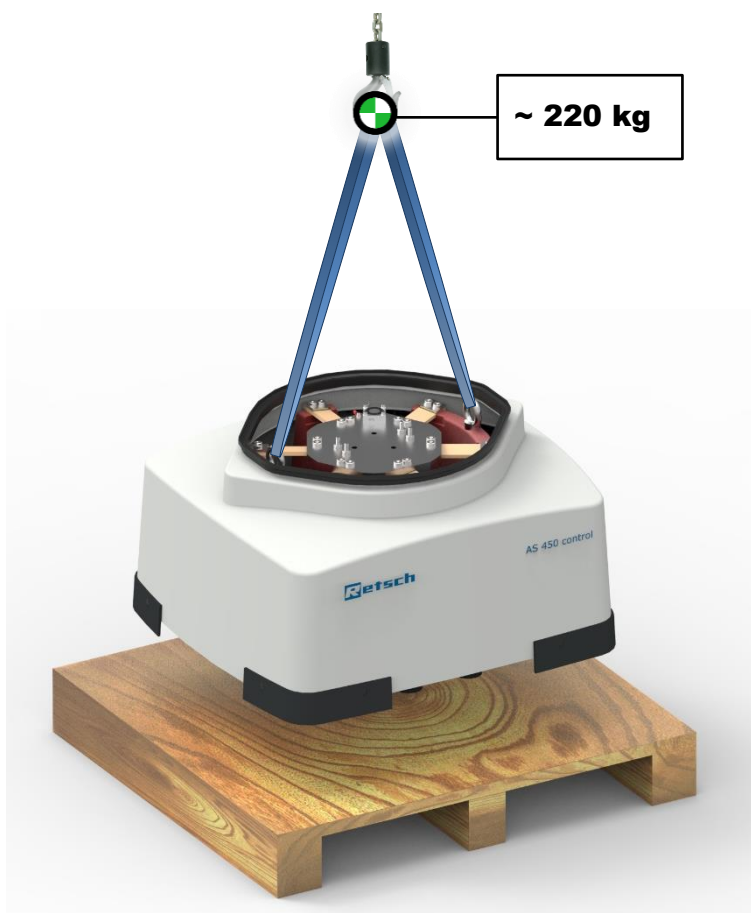



Fig. 2: Lifting the device

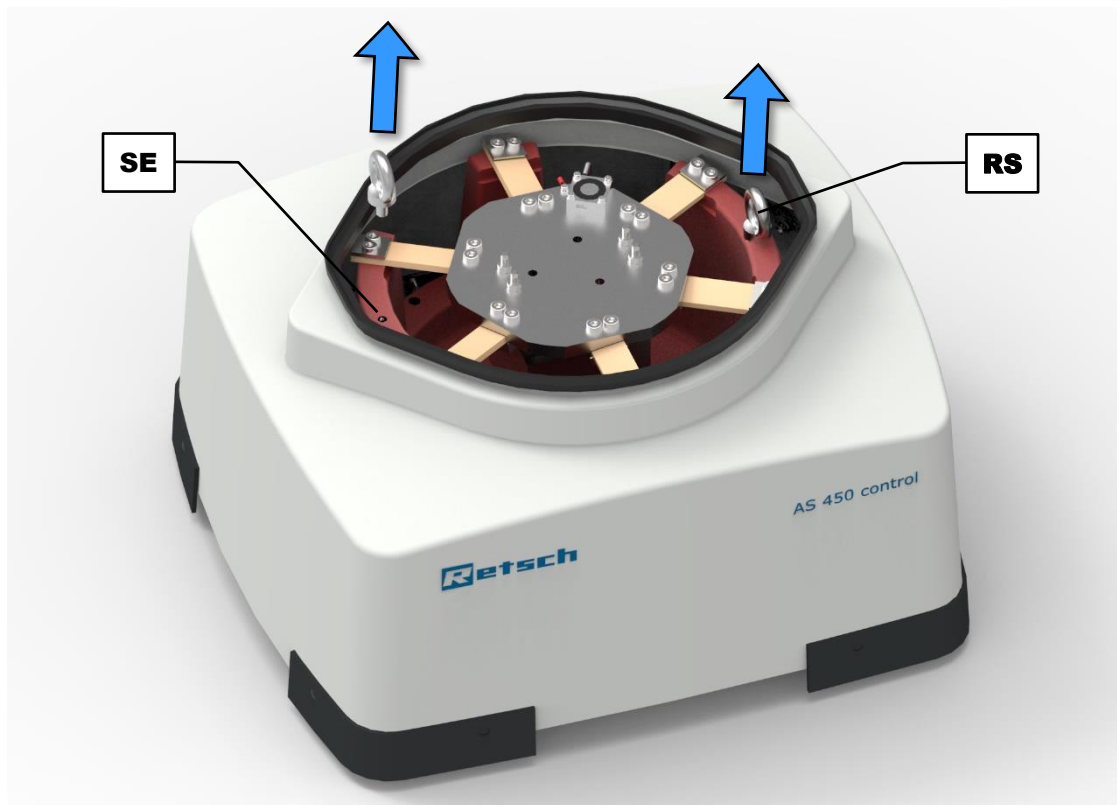


Fig. 3: Removing the transportation aid

The transportation aid consists of two eye bolts (**RS**).

- ⇒ After the installation of the device at the installation site, loosen the eye bolts (**RS**) on both sides of the oscillating unit (**SE**) and remove them.
- ⇒ Keep the transportation aid for later transport.

⚠ CAUTION The weight without sieve stack and sieve clamping unit amounts approx. 220 kg. The device may only be lifted with suitable hoist, that is designed for the weight.

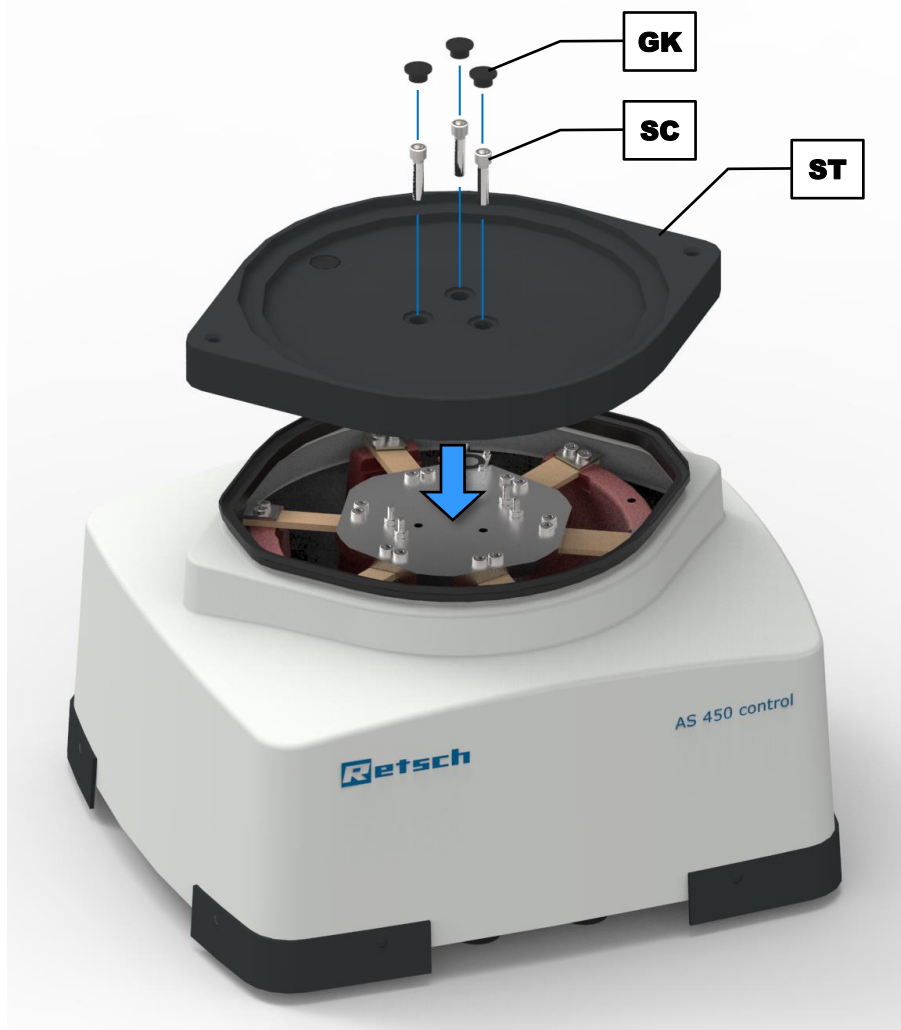


Fig. 4: Mounting the sieve plate

- ⇒ Place the sieve plate (**ST**) on the device. Ensure the correct alignment.
- ⇒ Tighten the sieve plate (**ST**) by means of the three supplied hexagon socket head screws (**SC**). The required torque is approx. 40 N·m.
- ⇒ Place the three supplied protective caps (**GK**) on the hexagon socket head screws and press them down firmly.


5 First Commissioning

⚠
WARNING

W4.0002

Danger to life through electric shock
Damaged power cable

- Operating the device with a damaged power cable or plug can lead to life-threatening injuries caused by an electric shock.
- **Before operating the device, check the power cable and plug for damage.**
- **Never operate the device with damaged power cable or plug!**



NOTICE

N10.0002

Setting up the device
Disconnecting the device from the mains

- A separation of the device from the mains must be possible at any time.
- **Set up the device in such a way, that the connection for the power cable is always easily accessible.**

NOTICE

N11.0004

Setting up the device
Vibrations during operation

- Depending on the operating mode of the device, slight vibrations may occur.
- **Set up the device only on a vibration-free, plane and stable surface.**

Before first commissioning the sieve clamping unit must be installed, and the control unit must be connected.

The AS 450 control is suitable for test sieves with an outer diameter of 400 mm and 450 mm. Up to 13 fractions (12 test sieves plus collecting pan with 400 mm outer diameter), or 9 fractions (8 test sieves plus collecting pan with 450 mm outer diameter) can be clamped.

NOTICE A high number of test sieves can significantly increase the total weight of the load (sieve stack and sample material). Make sure not to exceed the maximum payload of 75 kg.

Different sieve clamping units and lids are available for the test sieves.

The following sieve clamping units are available for test sieves with an outer diameter of 400 mm and 450 mm (18"):

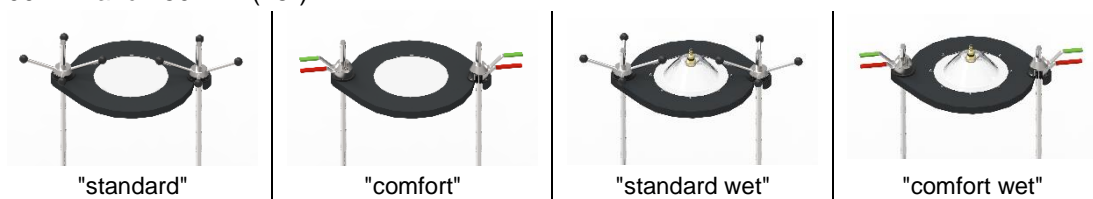


Fig. 5: Types of the sieve clamping unit

5.1 Sieve Clamping Unit "standard"

- ⇒ Screw both threaded rods (**A**) into the designated threaded holes (**SB**) in the sieve plate (**ST**).
- ⇒ Firmly tighten the threaded rods (**A**) by means of a 19 mm open-end wrench. The required torque is approx. 30 N·m.
- ⇒ Place the desired [sieve stack](#) including the sample material centrally on the sieve plate (**ST**).
- ⇒ Lay the clamping lid (**D**) over the threaded rods (**A**) onto the top test sieve. The clamping lid is orientated so that the peripheral edge surrounds the test sieves and the threaded rods (**A**) are up to the stop of the guide grooves (**D1**).
- ⇒ Screw the fixing nut (**B**) down the threaded rod (**A**) onto the clamping lid.
- ⇒ Tighten the fixing nuts hand-tight.

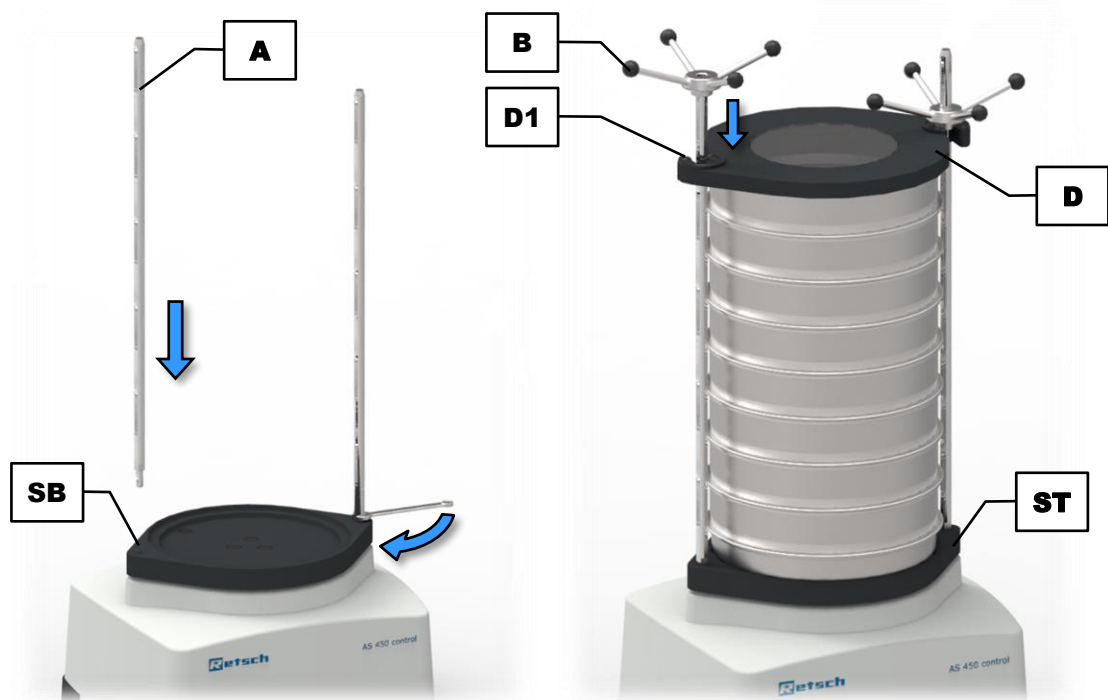


Fig. 6: Installation of the sieve clamping unit "standard"

5.2 Sieve Clamping Unit "comfort"

- ⇒ Screw both threaded rods (**A**) into the designated threaded holes (**SB**) in the sieve plate (**ST**).
- ⇒ Firmly tighten the threaded rods (**A**) by means of a 19 mm open-end wrench. The required torque is approx. 30 N·m.
- ⇒ Place the desired [sieve stack](#) including the sample material centrally on the sieve plate (**ST**).
- ⇒ Lay the clamping lid (**D**) over the threaded rods (**A**) onto the top test sieve. The clamping lid is orientated so that the peripheral edge surrounds the test sieves and the threaded rods (**A**) are up to the stop of the guide grooves (**D1**).

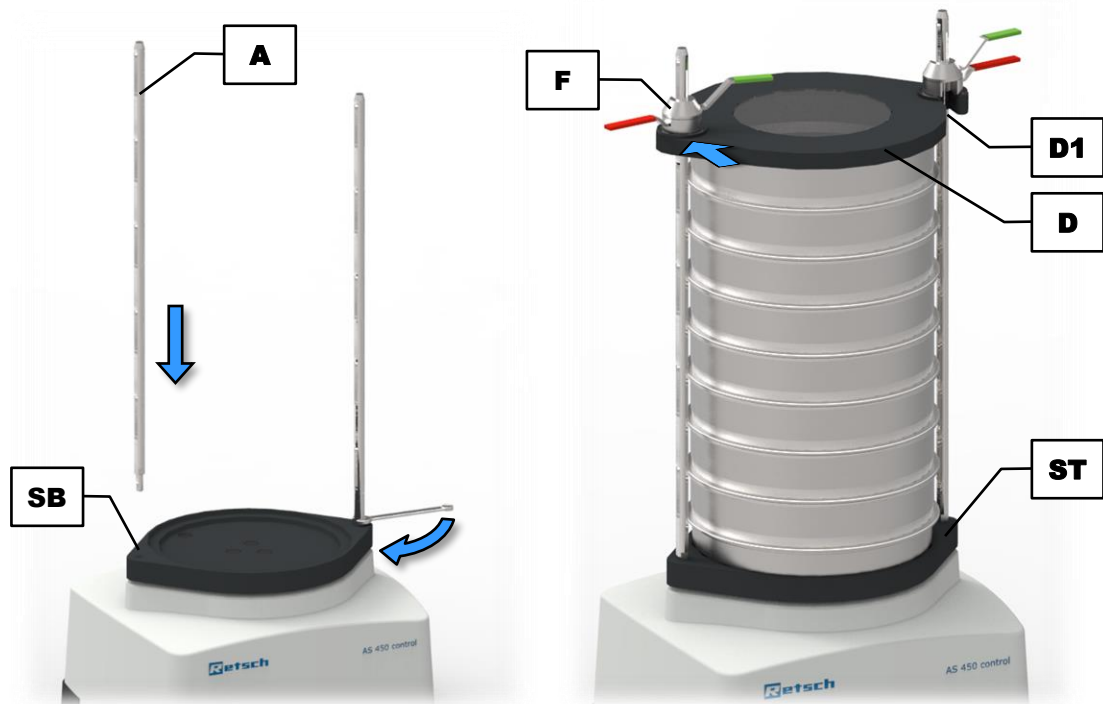
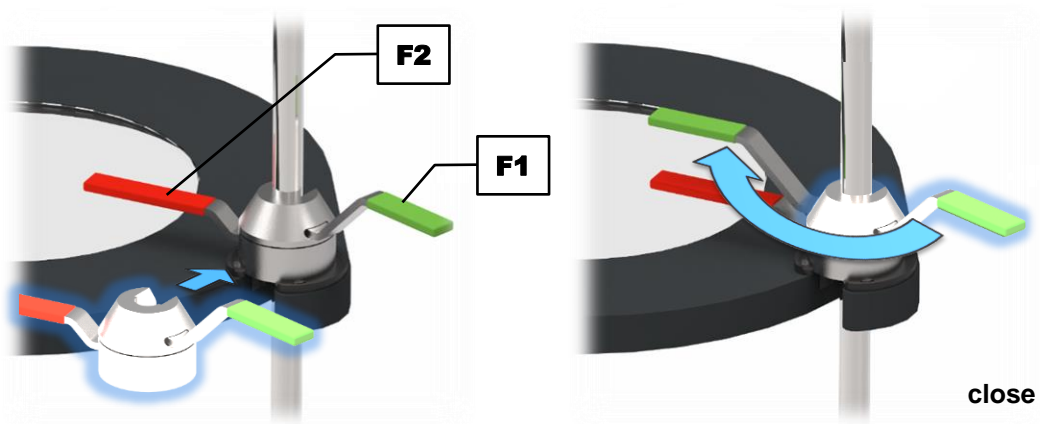


Fig. 7: Installation of the sieve clamping unit "comfort"

- ⇒ Open the quick clamping unit (F), i.e. the green (F1) and the red (F2) quick clamping levers are situated opposite of each other.
- ⇒ Place the opened quick clamping unit (F) on the threaded rod (A) directly above the clamping lid (D). Hereby, the threaded rod (A) must be enclosed by the quick clamping unit (F).



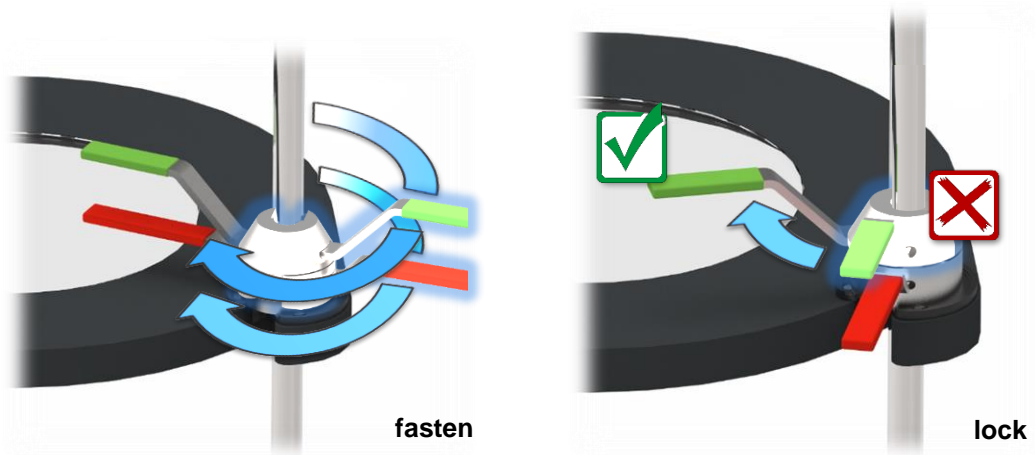


Fig. 8: Installation of the quick clamping unit

- ⇒ Close the quick clamping unit (F) by turning the green quick clamping lever (F1) until it is positioned above the red quick clamping lever (F2).
- ⇒ Turn both quick clamping levers (F1) and (F2) of the closed quick clamping unit (F) simultaneously clockwise, until the quick clamping unit (F) is hand-tightened.
- ⇒ Repeat these steps for the second quick clamping unit (F).
- ⇒ Fasten (lock) both quick clamping units (F) by turning the green quick clamping levers (F1) further clockwise as far as possible.

CAUTION

C6.0012

Contusions and bruises

Overturning of the sieve stack

- The sieve stack can overturn and cause personal injury.
- **Only operate the device with securely clamped sieve stack.**

5.3 Connecting the Control Unit

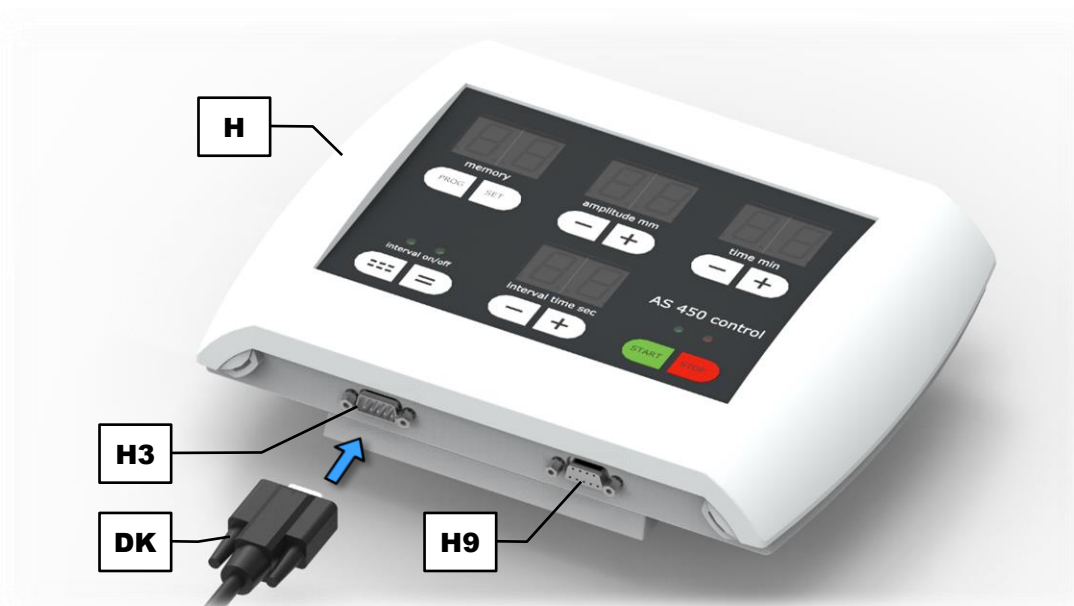


Fig. 9: Connecting the control unit

- ⇒ Connect the socket of the supplied data cable (**DK**) to the RS232 interface (**H3**) on the control unit (**H**).
- ⇒ Connect the plug of the supplied data cable (**DK**) to the RS232 interface (**O**) on the backside of the AS 450 control.
- ⇒ Secure both, the plug and the socket of the supplied data cable (**DK**) with the respective screws.

Via the RS232 interface (**H9**), a connection to a PC can be established. With the aid of the optionally available evaluation software EasySieve®, the AS 450 control can then be externally controlled and the results subsequently evaluated automatically (→ Chapter "[EasySieve®](#)").

5.4 Wall Mounting of the Control Unit

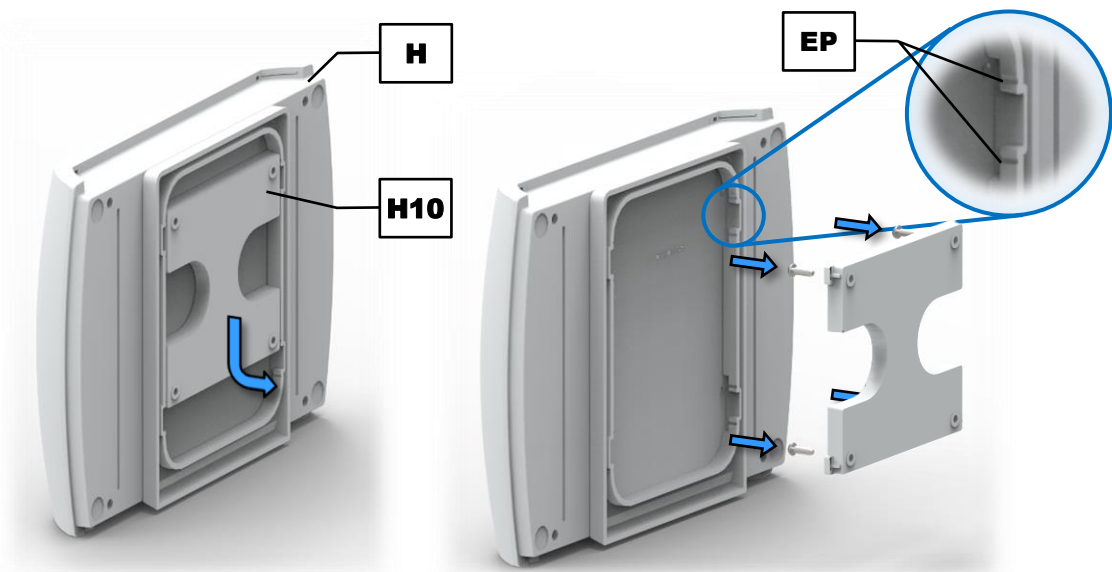


Fig. 10: Wall mounting of the control unit

The wall mounting plate (**H10**) is located on the backside of the control unit (**H**). The wall mounting plate has an upper and a lower latching position (**EP**) on the control unit.

- ⇒ Press the wall mounting plate (**H10**) upwards or downwards, depending on its latching position. The wall mounting plate can be removed from the control unit (**H**) when in mid position.

Four M5 screws with the appropriate dowels are required (not included in the scope of delivery) for the wall mounting.

- ⇒ Use the four screws to secure the wall mounting plate (**H10**) at the desired position on the wall. Make sure, that the flat surface of the wall mounting plate faces the wall.
- ⇒ Mount the control unit (**H**) on the wall by locking the wall mounting plate (**H10**) either in the upper or lower latching position (**EP**).

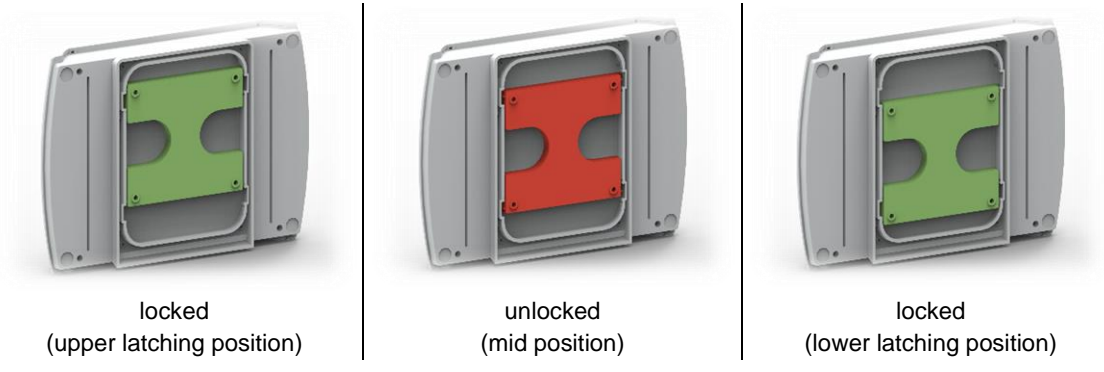


Fig. 11: Mounting the control unit

6 Operating the Device

6.1 Use of the Device for the Intended Purpose

⚠ CAUTION

C7.0005

Risk of injury

Potentially explosive atmosphere

- The device is not suitable for use in potentially explosive atmospheres. Operating the device in a potentially explosive atmosphere can lead to injuries caused by an explosion or fire.
- **Never operate the device in a potentially explosive atmosphere!**

⚠ CAUTION

C8.0006

Risk of injury

Sample material that is harmful to health

- Sample material that is harmful to health can injure people (illness, contamination).
- **Use suitable extraction systems with sample material that is harmful to health.**
- **Use suitable personal protective equipment with sample material that is harmful to health.**
- **Take note of the safety data sheets for the sample material.**



⚠ CAUTION

C9.0003

Risk of explosion or fire

Changing sample properties

- The properties and therefore also the hazardousness of the sample can alter during the sieving process.
- **Do not use any substances in this device which carry the risk of explosion or fire.**
- **Observe the material safety data sheets of the sample material.**



This Vibratory Sieve Shaker of the Retsch GmbH is a laboratory device. It is suitable for both, dry and wet sieving of free-flowing, disperse materials in the grain size range from 25 µm to 125 mm.

The particle size distribution of soils, building materials, chemicals, fertilizers, fillers, grains, coffee, plastics, flour, metal powders, minerals, nuts, seeds, sand, washing powder, cement clinker and many other substances can be easily and quickly analysed.


The Vibratory Sieve Shaker of the Retsch GmbH is successfully deployed in almost all areas of industry and research within the scope of quality control, especially where there are high demands regarding easy operability, speed, precision and reproducibility.

The AS 450 control is specially designed for test sieves with an outer diameter of 400 mm and 450 mm (18"). For an optimum measurement result it is recommended to exclusively use test sieves from Retsch GmbH.

⚠ WARNING W5.0010

Handling of food, pharmaceutical and cosmetic products
Analysed products

- Food, pharmaceutical and cosmetic products, which were analysed with the device must not be consumed, used or circulated.
- **Dispose these substances in accordance with the applicable regulations.**



NOTICE N12.0007

Range of application of the device
Long-term operation

- This laboratory device is designed for eight-hour single-shift operation with a duty cycle of 30 %.
- **This device may not be used as a production machine nor is it intended for continuous operation.**

6.2 Principle of Operation

The AS 450 control performs a vibratory sieving, where the sample material is thrown upwards by the vibrations of the sieve bottom and subsequently falls back down onto the sieve mesh fabric due to gravitation forces. Thereby, the sample material is subjected to a three-dimensional movement, i.e. a horizontal circular motion superimposes the vertical throwing motion. Hence, the sample material is spread uniformly across the entire surface of the sieve bottom, whereas the particles are subjected to an acceleration in vertical direction. In this process, they perform free rotations and are compared with the mesh sizes when falling back down statistically orientated. In the Vibratory Sieve Shaker of the Retsch GmbH, an electromagnetic drive sets a spring-mass system in motion and transfers the oscillations to the sieve stack. The amplitude can be adjusted within a few millimetres.

6.3 Views of the device

6.3.1 Front

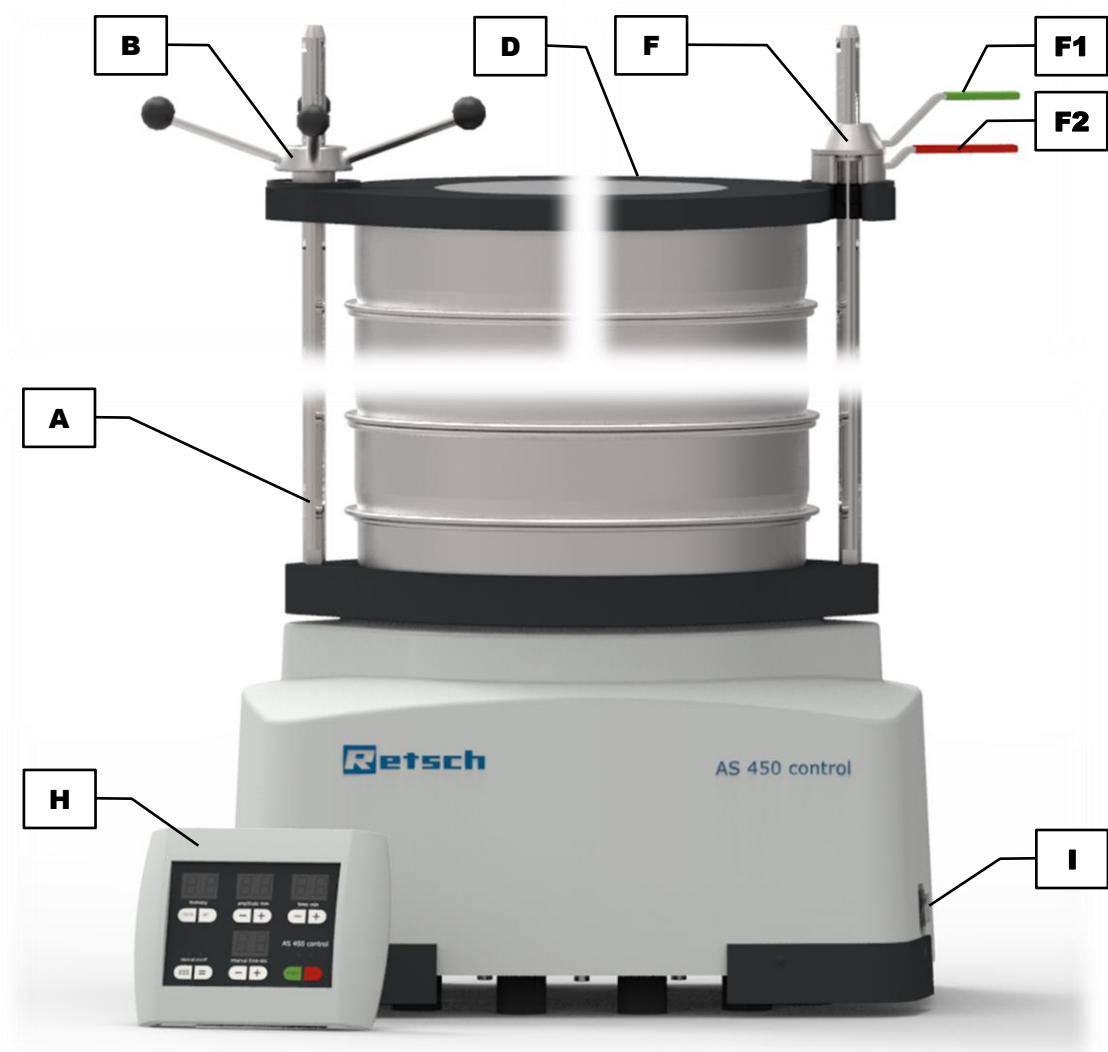


Fig. 12: Front view of the device with different sieve clamping units

Element	Description	Function
A	Threaded rod	Fixes the sieve stack together with the clamping lid (D) and the fixing nut (B) or the quick clamping unit (F)
B	Fixing nut "standard"	Fixes the sieve stack together with the clamping lid (D) and the threaded rod (A)
D	Clamping lid	Covers the top test sieve and fixes the sieve stack together with the fixing nut (B) or quick clamping unit (F) and the threaded rod (A)
F	Quick clamping unit "comfort"	Fixes the sieve stack together with the clamping lid (D) and the threaded rod (A)
F1	Quick clamping lever green	Clamps or releases the quick clamping unit (F) and thus, the sieve stack

Element	Description	Function
F2	Quick clamping lever red	Secures the quick clamping unit (F) and thus, the sieve stack
H	Control unit	Operation of the device
I	Mains switch	Switches the device on and off, disconnects the device from the mains

6.3.2 Back

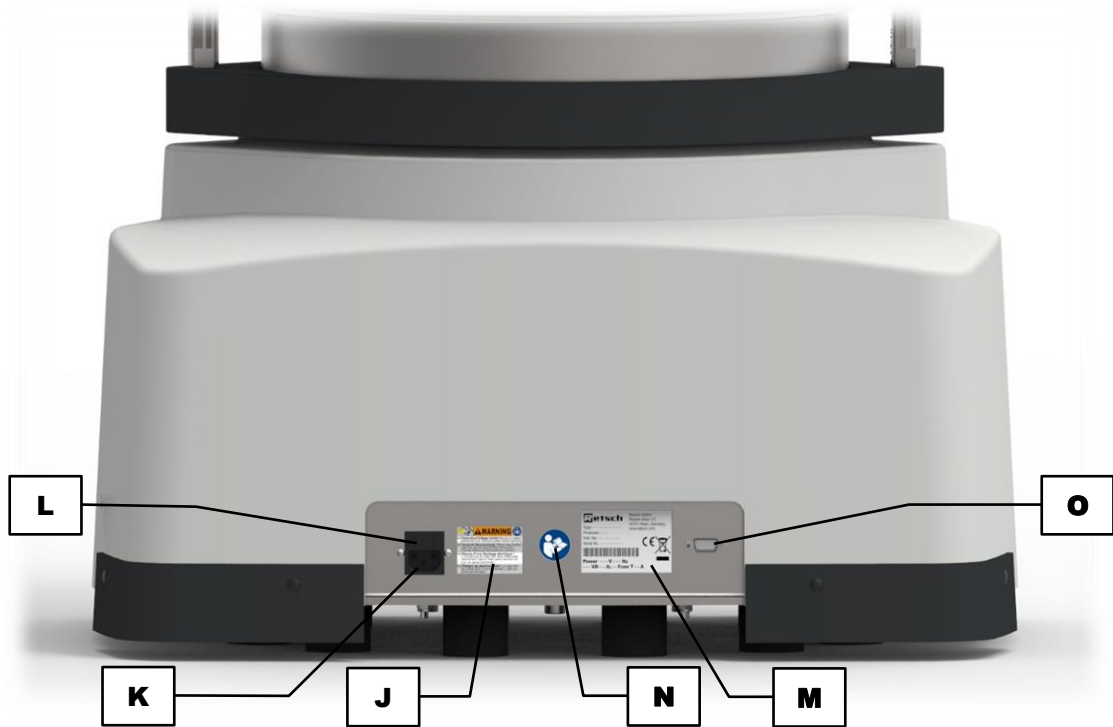


Fig. 13: Back view of the device

Element	Description	Function
J	Warning sign "Disconnect from the mains"	Warning of electric shock
K	Mains connection	Connection for the power cable
L	Fuse drawer	Contains the fuses protecting against overvoltage (fuse: 6.3 A delay-action at 220 – 240 V or 8 A delay-action at 100 – 120 V)
M	Type plate	Lists, among others, the voltage type, the serial number and the type of the device
N	Sticker "Manual"	Reminds to read the manual
O	RS232 interface	Connection for the control unit (H)

6.4 Switching On / Off




⇒ Turn on the AS 450 control with the mains switch (**I**) on the right side of the device.

When the device is switched off, it is completely disconnected from the mains.

Setting mode:

After switching on, the device is in the setting mode. The displays "time" (H5) and "amplitude" (H4) show the last used values. The display "memory" (H8) indicates "on" and the interval function is off.

Standby mode:

By pressing the  button (H1) after power on, the device can be put into standby mode. In this mode, only the LED of the  button (H1) is lit. The display "memory" (H8) indicates "on" and the interval function is off. All other displays are off. Except for the  button (H2), all buttons are inoperable.

6.5 Selection of the Test Sieves

The selection of the test sieves depends on the sample quantity as well as the particle size distribution. The gradation of mesh sizes and accordingly the measurement points should be selected in such a way that the complete particle size range of the sample is covered at regular intervals. The wider the particle size range, the more test sieves should be used.

NOTICE At least three test sieves and the collecting pan must be clamped (→ Chapter "[Amplitudes in Dependence on the Load](#)").

6.6 Performing a Sieving

- ⇒ Determine the empty weights of the test sieves and the collecting pan.
- ⇒ Place the sieve stack with **increasing** mesh size on the collecting pan.
- ① Each test sieve is provided with an O-ring, which serves as a seal to prevent dust emission during the sieving.
- ⇒ Weigh the sample and put it on the uppermost test sieve (biggest mesh size). Make sure not to exceed the [maximum feed quantity](#).
- ⇒ Place the complete sieve stack centrally on the device and clamp the sieve stack (→ Chapter "[Sieve Clamping Unit "standard"](#)" or "[Sieve Clamping Unit "comfort"](#)").
- ⇒ Set the optimum amplitude value and sieving time (→ Chapter "[Controlling the Device](#)").
- ⇒ Start the sieving process.
- ⇒ After the end of the sieving process, weigh the individual test sieves and the collecting pan including the particle size fractions present therein.
- ⇒ Determine the mass of the particle size fractions (weight after the sieving less the respective empty weight).
- ① The evaluation software "[EasySieve](#)[®]" automatically records the weights and allows for a quick and simple evaluation of the sieve analysis. When the device is controlled via EasySieve[®], "ES" is indicated in the display "memory" (H8). A detailed description can be found in the separate manual of the software.

7 Controlling the Device

7.1 Operating Controls, Displays and Functions

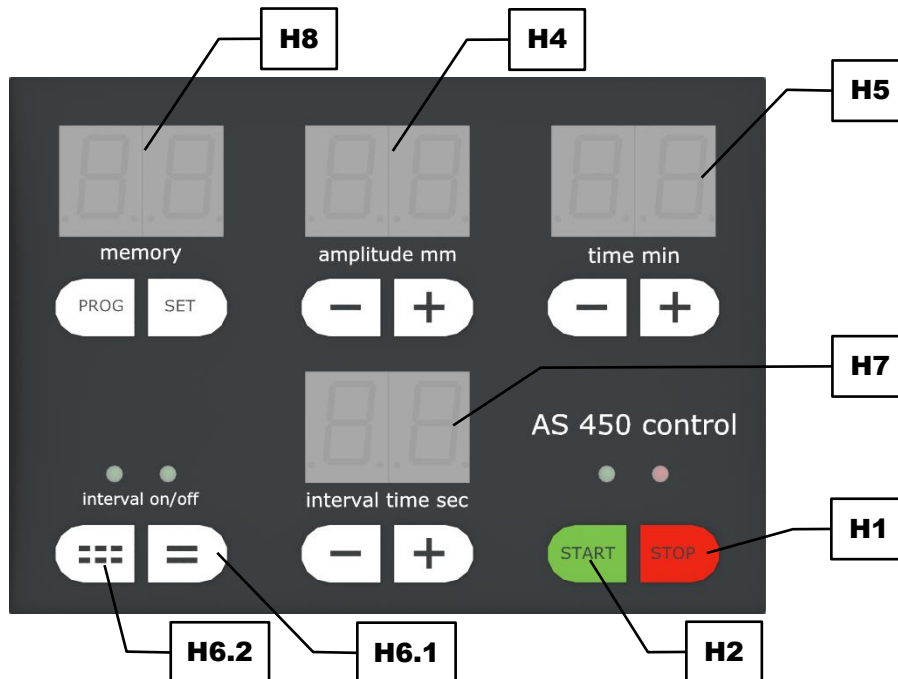



Fig. 14: Operating controls and functions

Element	Description	Function
H1	STOP	Stops the sieving process. In standby mode, the red LED is lit
H2	START	Starts the sieving process. During operation, the green LED is lit
H4	Amplitude setting	Decreases or increases the amplitude by pressing the "-" or "+" button, respectively in the range of 0.20 to 2.20 mm or 1 to 7.1 g
H5	Time setting	Reduces or extends the sieving time by pressing the "-" or "+" button, respectively in the range of 1 to 99 minutes
H6.1	Interval off	Switches the device to continuous operation. During the continuous operation, the green LED is lit
H6.2	Interval on	Switches the device to interval operation. During the interval operation, the green LED is lit
H7	Interval setting	Reduces or extends the sieving time between the interval pauses by pressing the "-" or "+" button, respectively in the range of 10 to 99 seconds
H8	Programme setting	Allows for the saving, editing and selection of up to 9 programmes

7.1.1 Start Process

⇒ To start the sieving process in the [setting mode](#), press the  button (**H2**).

- ⇒ If the device is in [standby mode](#), press the  button (**H2**) **twice** to start the sieving process.

The green LED lights up and the sieving process is started with the current settings.

7.1.2 Stop Process

The sieving process will stop automatically after the set process time has elapsed. However, the sieving process can be stopped manually at any time.

- ⇒ Press the  button (**H1**) **twice** to stop the sieving process.


7.2 Pause Process

The sieving process will stop automatically after the set process time has elapsed. However, the sieving process can be interrupted manually at any time.


- ⇒ Press the  button (**H1**) **once** to interrupt the sieving process.

The process time is stopped.

Continue the process:

- ⇒ Press the  button (**H2**) to continue with the sieving process.

End the process:

- ⇒ Press the  button (**H1**) to end the sieving process.

7.3 Amplitude

The amplitude display (**H4**) shows the set amplitude value depending on the setting either in mm or g (acceleration of gravity). The amplitude value is adjustable between 0.2 mm (1 g) and 2.2 mm (7.1 g) in steps of 0.1 mm (0.1 g). When the device is switched on, the last used amplitude is preset.

- ⇒ Press the "+" or "-" button to set the desired amplitude.
⇒ Press and hold the "-" or "+" button to increased or decreased the amplitude in fast steps.

The amplitude can also be changed during operation by pressing the "+" or "-" button. An exceeding or falling below of 2.2 mm or 0.2 mm respectively is not possible.

The amplitude can either be displayed as the height of lift in mm or as a multiple of gravity g (9.81 m/s²), the so called gravitational acceleration.

- ⇒ Simultaneously press the "+" and "-" button to toggle the display between "mm" and "g".

During the sieving process, the amplitude is kept constant within a predetermined tolerance of 0.1 mm.

NOTICE

N13.0008

Movement of the sieve stack at the beginning of the sieving process

Insufficient clamping of the sieve stack

- During the sieving process, the sieve stack can rotate and move over the sieve plate.
- The amplitude regulation cannot find a stable oscillating system, which leads to fluctuating amplitude values.
- **Observe the sieve stack at the beginning of the sieving process.**
- **If a movement of the sieve stack can be noticed, pause the sieving process and retighten the sieve stack again.**

7.3.1 Amplitudes in Dependence on the Load

The AS 450 control is a resonance sieving machine whose attainable amplitude is depending on the load. In this respect, the mass (sieve stack and sieve clamping unit) fixed to the sieve plate (**ST**) plays a primary role.

Only the amplitudes specified within the following load diagram can be achieved. The diagram is to be seen as a guideline for the voltage rated on the type plate (**M**). Mains voltage fluctuations or mains voltage deviations lead to increased tolerances.

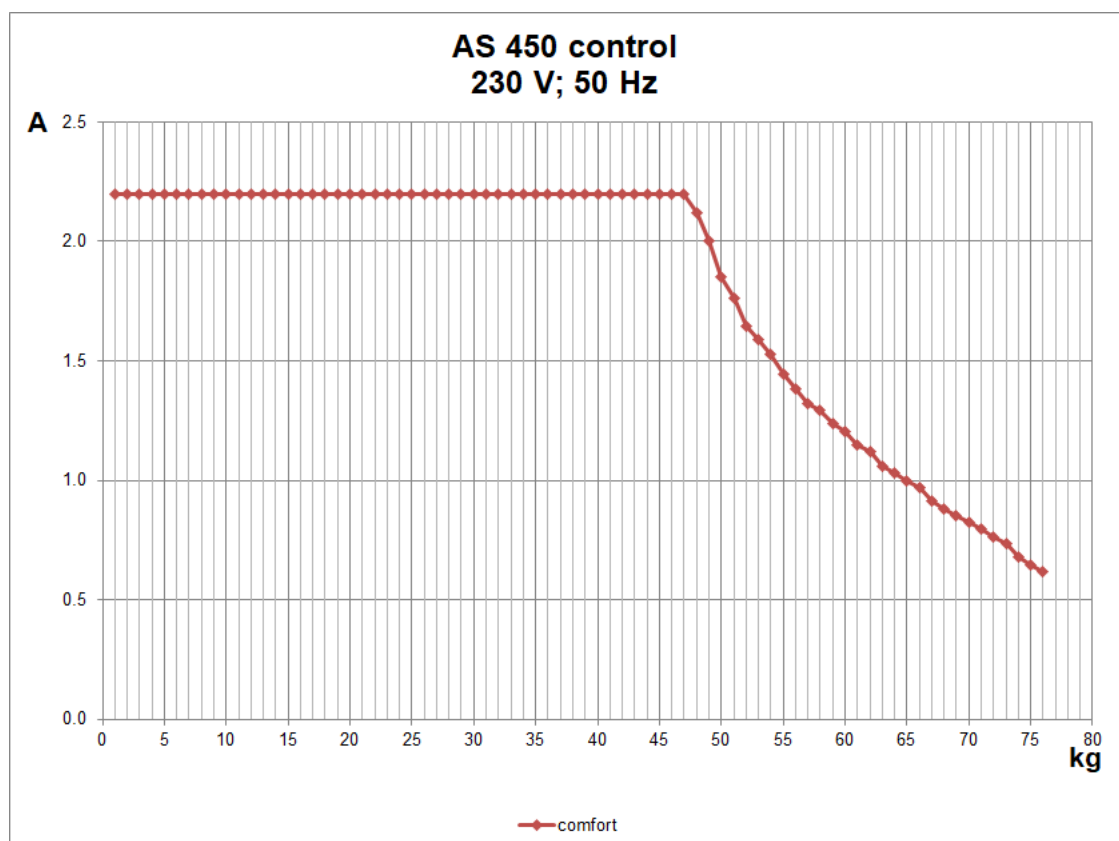


Fig. 15: Load diagram for the sieve clamping unit "comfort"

The diagram illustrates the amplitude "A" in millimetre in dependence of the load "kg" (sieve stack mass) in kilogramme. The tolerance of the sieve stack mass amounts $\pm 5\%$. The best results are generally achieved with an amplitude around 1.8 mm.

Example 1:

Type:	230 V; 50 Hz
Sieve stack mass:	20 kg
Sieve clamping unit:	"comfort"

Under these conditions the maximum attainable amplitude amounts 2.2 mm.

Example 2:

Type:	230 V; 50 Hz
Sieve stack mass:	50 kg
Sieve clamping unit:	"comfort"

Under these conditions the maximum attainable amplitude amounts ~ 1.9 mm.

7.4 Time

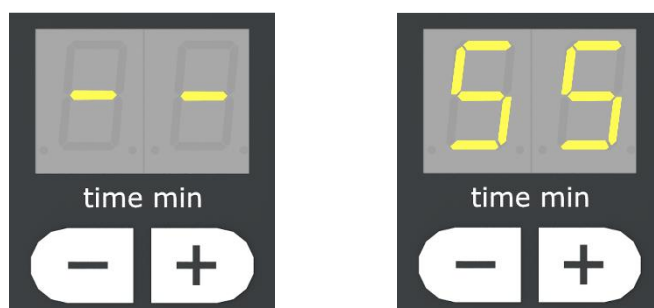


Fig. 16: Time setting for continuous operation (left) or with a process time (right)

The AS 450 control can be operated either in continuous operation or for a certain time between 1 and 99 minutes. When the device is switched on, the last used setting is displayed.

- ⇒ Press the "+" or "-" button of the time display (**H5**) to set the desired process time.
- ⇒ Press and hold the "+" or "-" button to extend or reduce the process time in fast steps.
- ⇒ To change to the continuous operation, fall below the duration of 1 min by pressing the "-" button, or exceed the duration of 99 min by pressing the "+" button. The time display (**H5**) now indicates "--".

The process time can also be changed during operation by pressing the "+" or "-" button.

7.5 Optimisation of Time and Amplitude

The settings of the optimum sieving time and amplitude depend on the sample material. These settings have a substantial influence on the measurement result. Generally, national and international standards, internal regulations and standards provide detailed information on product-specific sieve analyses and the associated sieving parameters. If such basic information cannot be obtained, the sieving time and amplitude must be determined experimentally.

With the AS 450 control the amplitude is defined as the total lifting height (**SH**) of the test sieve. For example, with a set amplitude of 1.2 mm, the test sieve is displaced in the range of -0.6 mm and +0.6 mm around the zero point (= stationary sieve plate (**ST**)).

An **optimum amplitude** has been found, when a state of statistical resonance is being reached during the sieving process. Then, the particles have the biggest probability of passing, as the

throw time of a particle corresponds to the oscillation period of the test sieve. In this case, the particle (**PA1**) will be moved with a different orientation to a different mesh every time the test sieve (**SH**) lifts. At too low amplitudes, the particles (**PA2**) do not lift off high enough from the sieve mesh fabric, and are therefore not able to orientate freely and move freely over the sieve mesh fabric. At too high amplitudes, the particles (**PA3**) are thrown up very high, and thus have fewer opportunities to compare themselves with the sieve meshes. The best results are generally achieved with an amplitude around 1.8 mm.

The **optimum sieving time** is in accordance with DIN 66165 achieved, if less than 0.1 % of the feed quantity passes the test sieve after one minute of sieving duration. In practice, the individual test sieves are weighed after the sieving process including the respective particle size fraction. Then, the sieve stack is sieved again for one minute. The weights of the individual test sieves of the second weighing must not differ substantially from those of the first weighing.

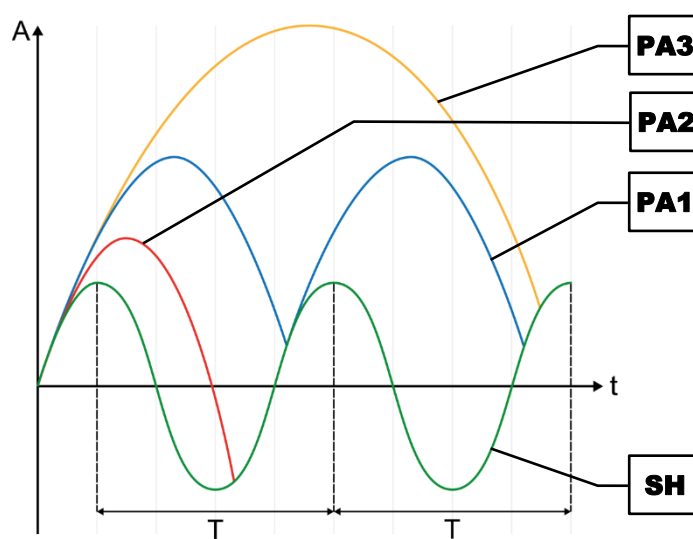
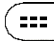



Fig. 17: Movement of the particles on the test sieve

7.6 Interval

- ⇒ Press the  button (**H6.2**) to change to the interval operation. The corresponding green LED is lit and the display "interval time" (**H7**) is enabled.
- ⇒ Press the  button (**H6.1**) again, to change back to the continuous operation. The corresponding green LED is lit and the display "interval time" (**H7**) is disabled.

During the interval operation, the sieving process is periodically interrupted for approx. one second. In interval operation, the interval times (pause times) are included in the displayed process time (**H5**). The interval operation cannot be switched on or off during the sieving process.

7.6.1 Interval Time

With the interval time (**H7**) the time of the sieving process between the interval pauses can be freely selected between 10 and 99 seconds. The pause time of approx. one second is not changeable. The display of the interval settings is only active with the interval function turned on.

- ⇒ Press the "+" or "-" button of the interval time (**H7**) to set the desired sieving time between the interval pauses.
- ⇒ Press and hold the "+" or "-" button to extend or reduce the interval time in fast steps.

On exceeding 99 seconds, the display changes back to 10 seconds. When going below 10 seconds, the display changes to 99 seconds. The interval time cannot be changed during the sieving process.

7.7 Programme Mode

The AS 450 control allows for the saving and recalling of up to 9 parameter sets. The programme settings can only be edited in the [setting mode](#).

The currently selected programme is displayed in the display "memory" (H8). If the display indicates "on", no programme is selected and the device is in the manual mode.

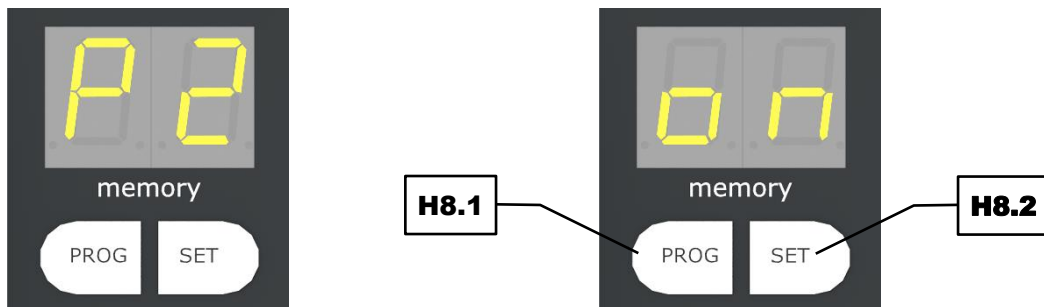











Fig. 18: Programme mode (left), manual mode (right)


7.7.1 Select a Program

- ⇒ Press the  button (H8.1) to navigate through the programme memory slots in ascending order and to select the desired programme.
- ⇒ Press the  button (H2) to start the sieving process in programme mode. All buttons, except for the  button (H1) are now locked.


Following the programme memory slot P9, "on" is again displayed in the display "memory" (H8) and the device is in manual mode. When a programme is selected, all buttons except for the  (H8.1),  (H8.2),  (H2) and  (H1) button are locked.

7.7.2 Edit a Program

- ⇒ Press the  button (H8.1) until the desired programme memory slot is displayed.
- ⇒ Press the  button (H8.2). All displays are now flashing.
- ⇒ Set the desired sieving parameters (amplitude, time, interval).

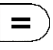

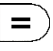

The programming can be cancelled by pressing the  button (H8.1). All settings are discarded.

7.7.3 Save a Programme


- ⇒ Press the  button (H8.2) to save the set sieving parameters in the selected programme memory slot. The displays stop flashing.

7.8 Signal Tone

The end of the sieving process is announced with an acoustic signal.

- ⇒ Simultaneously press the  (H6.1) and  (H1) button to turn off the signal tone. The process is confirmed by a signal tone.
- ⇒ Simultaneously press the  (H6.1) and  (H2) button to turn on the signal tone. The process is confirmed by a signal tone.

7.8.1 Operating Hours

- ⇒ Simultaneously press the  button (H8.2) and the "+" button of the interval time display (H7).

The display "memory" (H8) indicates "bS" (Betriebsstunden = operating hours). The complete runtime (corresponds to the accumulated sieving duration) of the device is displayed in the following three displays in hhhh:mm format:

- the minutes in mm format in the interval time display (H7)
- the hours in hhhh format split into the two displays "amplitude" (H4) and "time" (H5)

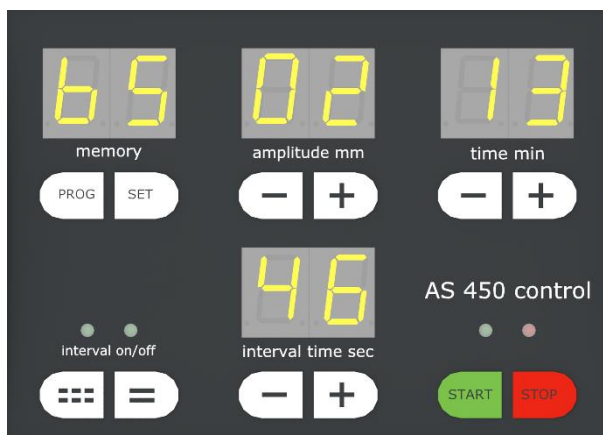




Fig. 19: Example of the display of the operating hours with a complete runtime of 213 hours and 46 minutes

All buttons, except for the  button (H1) are now locked.

- ⇒ Press the  button (H1) to exit the display of the operating hours.

7.8.2 Software Version

- ⇒ Simultaneously press the  button (H8.2) and the "-" button of the interval time display (H7).

The display "memory" (H8) indicates "S" (software). The current software version is displayed in the two displays "amplitude" (H4) and "time" (H5).

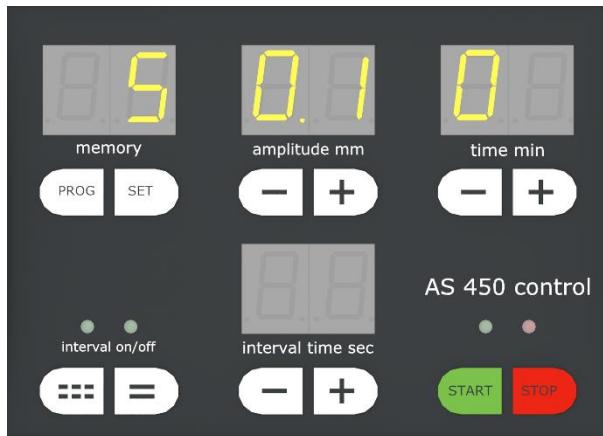



Fig. 20: Example of the display for the software version 0.10

All buttons, except for the  button (**H1**) are now locked.


⇒ Press the  button (**H1**) to exit the display of the software version.

8 Wet Sieving

⚠ WARNING W6.001

Danger to life through electric shock
Wet sieving


- An electric shock can cause burns, cardiac arrhythmia, respiratory arrest, as well as cardiac arrest.
- **Never operate the device in a water drain basin!**
- **Do not touch the device, if water has entered the interior!**
- **Always operate the device with a mains socket protected by a residual current circuit breaker (RCCB).**



⚠ WARNING W7.006

Danger to life through electric shock
Ingress of water into the control unit


- The control unit is not protected against the ingress of water. Water inside the control unit can cause a short circuit and an electric shock.
- **Make sure that the control unit does not come into contact with water!**
- **Do not touch the control unit, if water has entered the interior!**
- **Pull out the power plug!**



⚠ WARNING W8.008

Risk of death caused by an electric shock
Penetration of water if the mains plug is not fully plugged in

- If the IEC connector is not fully plugged in to the IEC appliance socket, water may enter the socket causing an electric shock.
- **Only operate the device with the IEC connector fully plugged in.**



NOTICE N14.0049

Damage to the sieve mesh fabric
Fluid retention during wet sieving

- Fluid retention can lead to overload and therefore to the damage or destruction of the sieve mesh fabric.
- **Observe the recommended flow rate.**
- **Always dose the quantity of liquid applied in such a way that no fluid retention can occur.**
- **Use venting rings, if necessary.**

8.1 Performing the Wet Sieving

Usually, sieving processes are carried out dry. However, when agglomerates, electrostatic charges or a high degree of fines impede the sieving process, either [sieving aids](#) can be used, or a wet sieving can be performed.

For wet sieving, a liquid, preferably water, is supplied to the sample material during the sieving process. A condition for wet sieving, however, is that the material to be sieved does not swell,

dissolve or otherwise change in the liquid. Wet sieving is particularly suitable for materials which are already in suspension and may not be dried.

In addition to the test sieves, a collecting pan (AB1) with an outlet (AB2) and a wet sieving lid (ND1) with spray nozzle (ND2) are required for wet sieving. During the sieving process liquid is introduced via the spray nozzle (ND2) situated on top of the upper most test sieve into the sieve stack and, subsequently, leaving it again together with the last fraction via the outlet (AB2) of the collecting pan (AB1).

- ⇒ Position the device in the vicinity of the drain point (e.g. drain in the floor). The distance between the outlet (AB2) and the drain point should not be too large.
- ⇒ Connect the spray nozzle (ND2) of the wet sieving lid (ND1) with the liquid supply (e.g. water tap). The inner diameter of the hose must be 13 mm.
- ⇒ Connect the outlet (AB2) of the collecting pan (AB1) with the drain point or a corresponding receptacle. The inner diameter of the hose must be 20 mm. Make sure that the drain point or the receptacle are located **below** the collecting pan (AB1) and that the hose has a **continuous slope down**.

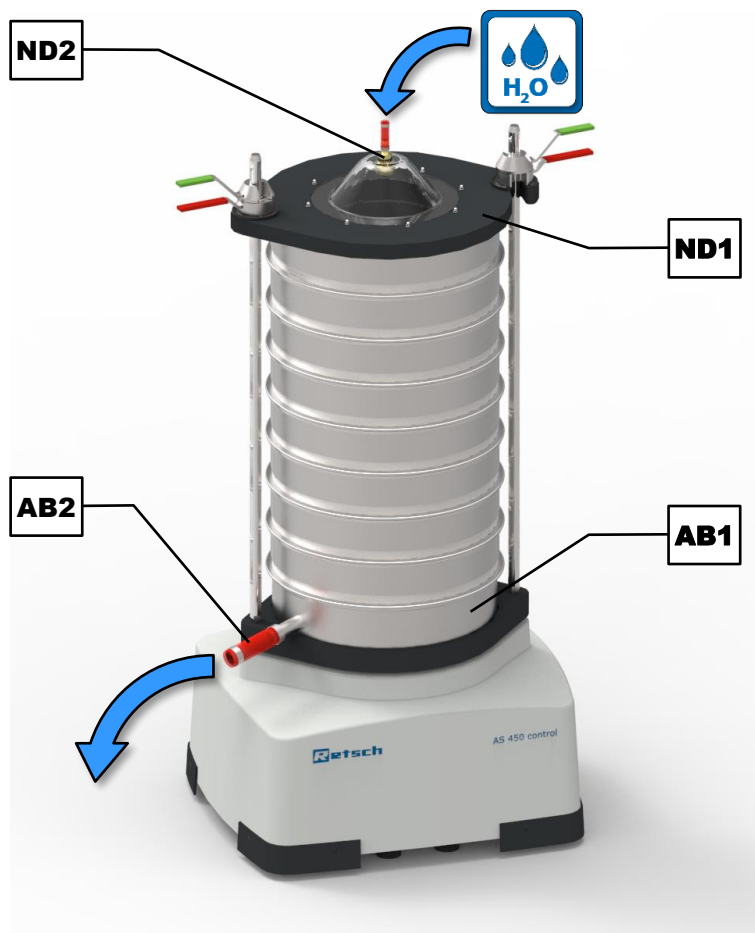


Fig. 21: Wet sieving

- ⇒ Suspend the sample material in a beaker containing the liquid intended for the wet sieving. To reduce the surface tension and to facilitate the screenings of the material later on, a few drops of surfactant may be added.
- ⇒ Moisten each test sieve with the liquid intended for wet sieving.
- ⇒ Place the sieve stack with **increasing** mesh size on the collecting pan with outlet.
- ⇒ Place venting rings (ER) between test sieves of mesh size < 100 µm to avoid air cushions.
- ⇒ Place the complete sieve stack centrally on the device.

- ⇒ Enter the sample suspension on the uppermost test sieve with the clamping lid open.
- ⇒ Clamp the sieve stack (→ Chapter "[Sieve Clamping Unit "comfort"](#)" or "[Sieve Clamping Unit "standard"](#)").
- ⇒ Set the optimum amplitude value and sieving time (→ [recommended parameters](#)).
- ⇒ Start the sieving process.
- ⇒ Turn on the liquid supply. The quantity of liquid applied should only be of such amount that the sieve mesh area is completely sprayed. A flow rate of 200 to 300 ml per sieve surface in dm² and minute is recommended (e.g. 3.2 to 4.8 litre per minute for a sieve diameter of 450 mm).
- ⇒ The sieving process is considered as terminated when the exiting liquid shows no turbidity anymore.



Fig. 22: Venting ring

If the smallest fraction, that leaves the collecting pan should also be weighted, it must be appropriately collected. After the sieving process, the individual fractions are transferred on suitable tared filters (paper filter) and dried in an oven at 80 °C until the weight remains constant.

NOTICE Used test sieves must be cleaned immediately after the sieving process (→ Chapter "[Cleaning of Test Sieves](#)"). Depending on the sample material flash rust can form in the sieve mesh fabric.

- ① The [load diagrams](#) are invalid for the wet sieving. Due to the non-defined quantity of liquid in the sieve stack, binding statements are not possible for the wet sieving.

Recommended parameters for wet sieving:

- Amplitude: 1.2 mm to 1.5 mm
- Interval operation: yes
- Time: 5 min

9 EasySieve®

EasySieve® is a software for particle size analysis and simplifies the manual evaluation in many respects. The software is able to automatically perform the required measuring and weighing processes – from determining the weights of the test sieves to evaluating the data.

The software is structured in a self-explanatory way and follows the logical chain of events involved in a particle size analysis. Therefore, the training period will be very short. The abundance of evaluation options provides absolute flexibility in adapting to demanding, individual tasks.

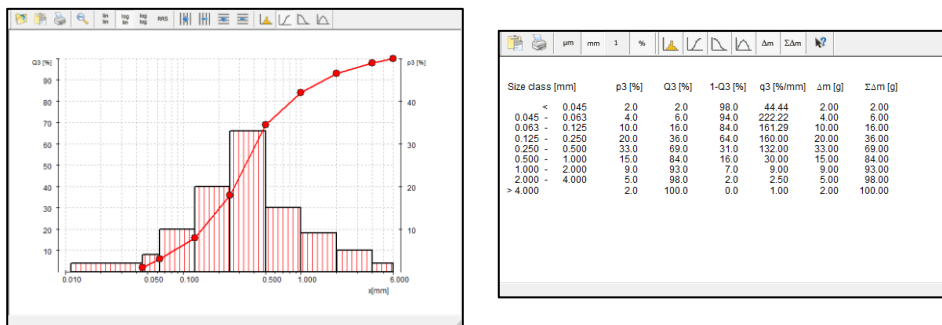


Fig. 23: Graphic and tabular presentation of the particle size analysis with EasySieve®

The software communicates with the scale and the AS 450 control and guides the user through the respective steps. Available parameters, as well as the characteristics to be calculated can be entered in various edit boxes. Routine parameters can be edited, saved and recalled at any time.

If a scale is connected, the corresponding data (empty weights of test sieves, back weights of loaded test sieves) can be transferred directly to EasySieve®. If no scale is connected, the input can also be entered manually.

The software calculates all standard particle distributions, as well as the representative characteristics of the particle size, and allows for the tabular and graphical presentation of results in a measurement report conforming to standards. Furthermore, the data can be exported to other software products (e.g. Microsoft Excel).

EasySieve® is also available as AuditTrail enabled version in compliance with 21CFR Part 11.

① In case the sieve shaker is controlled via the software "EasySieve CFR", warning messages may appear in the Audit Trail stating that no serial number, no software version and no calibration date is received. Moreover, it may indicate, that no error number can be received. This is correct because the sieving machine does not support the provision of this data. Therefore, these warning messages in the Audit Trail do not constitute a reason for corrective measures.

① A detailed description can be found in the separate manual of the software.

10 Error Messages and Information Notes



10.1 Error Messages

Error messages inform the user about detected device or programme errors. In the event of an error message, a fault has occurred, in which the operation of the device or the programme is automatically interrupted. Such faults must be resolved before next startup.

Error code	Description	Measures
E00	No errors in memory	⇒ The test of the memory has been completed without error.
E10	Drive overload	⇒ Switch off the main switch and wait for 30 s before switching on again. ⇒ If the error persists, contact service.
E20	Failure main board	⇒ Switch off the main switch and wait for 30 s before switching on again. ⇒ If the error persists, contact service.
E26	Failure frequency converter	⇒ Switch off the main switch and wait for 30 s before switching on again. ⇒ If the error persists, contact service.

10.2 Information Notes

Notices inform the user on specific device or programme processes. The operation of the device or programme may be interrupted briefly, but there is no fault. The information notice must be acknowledged by the user to continue the process. Information notices provide additional information for the user as an aid, but do not represent any device or programme errors.

Notice code	Description	Measures
bS	Display of the complete runtime in hhhh:mm	⇒ Press the  button to exit the display.
S	Display of the software version	⇒ Press the  button to exit the display.
ES	External control by EasySieve®	The device is controlled by the EasySieve® software installed on a PC. ⇒ Close the software to restore the manual control.

11 Return for Service and Maintenance



Fig. 1: Return form


The acceptance of devices and accessories of the Retsch GmbH for repair, maintenance or calibration can only be effected, if the return form including the decontamination declaration service has been correctly and fully completed.

- ⇒ Download the return form located in the download section "Miscellaneous" on the Retsch GmbH homepage (<http://www.retsch.com/downloads/miscellaneous/>).
- ⇒ When returning a device, attach the return form to the outside of the packaging.

In order to eliminate any health risk to the service technicians, Retsch GmbH reserves the right to refuse the acceptance and to return the respective delivery at the expense of the sender.

12 Cleaning, Wear and Maintenance



12.1 Cleaning


WARNING

W9.0003

Risk to life caused by an electric shock
Cleaning live parts with water

- Cleaning the device with water can lead to life-threatening injuries caused by an electric shock if the device has not been disconnected from the power supply.
- **Only carry out cleaning work on the device when it has been disconnected from the power supply.**
- **Use a cloth moistened with water for cleaning.**
- **Do not clean the device under running water!**

NOTICE

N15.0009

Damage to the housing and device
Use of organic solvents

- Organic solvents may damage plastic parts and the coating.
- **The use of organic solvents is not permitted.**

⇒ Clean the housing of the device with a damp cloth and if necessary, with a household cleaning agent. Pay attention that no water or cleaning agent enters the interior of the device.

12.1.1 Cleaning of Test Sieves

Test sieves are measuring instruments and should be treated with due care before, during and after the sieving process. It is recommended to clean new test sieves before the first use from possible preservative residues with ethanol or isopropanol and to store them in a dry, dust-free place when unused.

Before cleaning or drying the test sieves, the O-rings have to be removed. Before using and after the cleaning the test sieves should be visibly inspected for possible damages and impurities.

Near-mesh or clumped particles can be often removed dry after the sieving process by slightly tapping the test sieve upside down with the sieve frame on a table. For test sieves with mesh sizes > 500 µm a fine hair brush can be used to sweep over the outer side of the mesh fabric.

12.1.1.1 Cleaning of Test Sieves with Mesh Sizes > 500 µm

Coarse mesh fabrics with mesh sizes > 500 µm can be cleaned dry or wet easily and effectively with a hand brush with plastic bristles (at not too high applied pressure).

12.1.1.2 Cleaning of Test Sieves with Mesh Sizes < 500 µm

Test sieves with mesh sizes < 500 µm should generally only be cleaned in an ultrasonic cleaning-bath. As cleaning agent, water together with a standard surfactant is recommended. The cleaning in the ultrasonic bath usually takes two to three minutes. After that the test sieves

are thoroughly rinsed with water and dried. The cleaning with strong bases or acids is generally not recommended.

12.1.1.3 Drying of Test Sieves

Drying ovens of various sizes can be used for drying test sieves (drying temperature < 80 °C).

Additional information concerning ultrasonic cleaning-baths and drying ovens can be found on the Retsch GmbH homepage (<https://www.retsch.com>). Also ask for the free expert guide *Sieve Analysis – Taking a close look at quality*.

NOTICE

N16.0028

Damage of the sieve mesh fabric

Drying temperature > 80 °C

- At higher temperatures, especially fine metal wire meshes can become warped, leading to a reduced tension of the mesh fabric inside the sieve frame and hence, makes the test sieve less efficient during the sieving process.
- **The drying temperature for test sieves must not exceed 80 °C!**

12.2 Wear

Even with the proper handling of the test sieves, a wearing of the sieve mesh fabric depending on the frequency of the sieving operation and on the sample material is unavoidable. The test sieves should be regularly checked for wear and damage and be replaced if necessary.

Likewise, all existing sealing gaskets should be checked for wear on a regular basis and replaced if necessary.

⚠ CAUTION

C10.0013

Risk of injury

Improper repairs

- Unauthorised and improper repairs can cause injuries.
- **Repairs to the device may only be carried out by the Retsch GmbH , an authorised representative or by qualified service technicians.**
- **Do not carry out any unauthorised or improper repairs to the device!**


13 Maintenance

The AS 450 control is largely maintenance-free.

If wet sieving is executed, a quarterly examination for tightness of the fluid hoses should be performed.

If the AS 450 control is used in quality control, it should be regularly calibrated in accordance with DIN EN ISO 9000 ff. For this purpose please contact your local distributor or get in touch with Retsch GmbH directly.



13.1.1 Replacing the Fuses


WARNING

W10.0014

Risk to life caused by an electric shock
Exposed contacts

- Replacing the fuses without pulling out the mains plug can lead to life threatening injuries caused by an electric shock on contact with the fuse holder or the live contacts on the fuse.
- **Pull out the mains plug before replacing the fuses.**

NOTICE Depending on the mains supply different fuses are used. The correct electrical protection is listed on the type plate (**M**).

Voltage	Fuse
100 – 120 V	8 A delay-action
200 – 240 V	6.3 A delay-action

Two fuses are located in the fuse drawers (**L**) on the backside of the device. Fuses can be replaced by trained qualified personnel.

- ⇒ Unscrew the fuse drawers by means of a flat-bladed screwdriver.
- ⇒ Replace the defective fuses in the fuse drawers.
- ⇒ Screw the fuse drawers back in.

14 Accessories

Information on available accessories as well as the respective manuals are accessible directly on the Retsch GmbH homepage (<https://www.retsch.com>) under the heading "Downloads" of the device.

Information on wear parts and small accessories can be found in the Retsch GmbH general catalogue also available on the homepage.

In case of any questions concerning spare parts please contact the Retsch GmbH representative in your country, or Retsch GmbH directly.

14.1 Test Sieves

Decisive for the accuracy and reliability of the measurement result is, in addition to the reproducible operating Vibratory Sieve Shaker the quality of the test sieve. Test sieves of Retsch GmbH are high quality measuring instruments for which only mesh fabrics and perforated sheets of the corresponding standards are used. Each test sieve is tested five times and is given a serial number, as well as a quality certificate after the final check.



Fig. 24: Test sieves

The different versions of the test sieves of Retsch GmbH are supplied in accordance with all current national and international standards:

- available standards: DIN, ISO, ASTM, BS, NF, CGSB
- available diameters: 100 mm / 150 mm / 200 mm / 203 mm (8") / 305 mm (12") / 400 mm / 450 mm (18")
- available sieve surfaces: sieve mesh fabric (20 µm to 125 mm) and perforated screens (round, elongated or square holes) of stainless steel
- on request with an individual test certificate for the inspection of measuring and testing equipment monitoring according to ISO 9000 ff.

Among the various test sieves matching collecting pans, collecting pans with outlet, intermediate pans, intermediate rings, venting rings and sieve lids are available.

14.1.1 Certificate

Before delivery, each test sieve is optically surveyed according to the standards DIN ISO 3310-1 and ASTM E 11, and provided a certificate of compliance with the order.

On request, an additional acceptance test certificate with a calibration protocol can be provided, documenting the measurement results in tabular and graphical form, hence representing a calibration certificate with more detailed statistics.

14.1.2 Calibration Service

As a special service Retsch GmbH offers the calibration of the test sieves. All relevant information are recorded during the standard measuring process of the test sieve and confirmed in the required certificate.

14.2 Sieving Aids

NOTICE

N17.0027

Damage of the sieve mesh fabric

Use of mechanical sieving aids

- When using mechanical sieving aids, there is a danger that fine sieve mesh fabrics might be damaged.
- **Ensure that no overstretching of the sieve mesh fabric occurs due to overloading with sieving aids.**
- **If in doubt, please contact your local distributor or Retsch GmbH directly.**

By electrostatic and Van-der-Waals forces, as well as by fluid bridges, single particles can combine to form agglomerates. Since in this case not the individual primary particles, but particle collectives are measured, there is a distortion of the particle size distribution (a higher coarse fraction results). In order to prevent the formation of agglomerates or dissolve them, sieving aids can be used.

Mechanical sieving aids:

Mechanical sieving aids cause a destruction of agglomerates and dislodge wedged particles from the sieve meshes. Depending on the mesh size of the test sieve and the preselected amplitude, balls of agate, rubber, steatite or cubes of polyester urethane rubber, and nylon brushes can be used for this purpose.

NOTICE For very soft sample material, an undesired crushing of primary particles might occur.

Solid additives:

Solid additives, such as talcum or Aerosil® can be admixed to fatty, moist, sticky or oily sample materials. They attach themselves to the particle surface and counteract the formation of agglomerates. Their particle size is so small that they have no sustainable influence to the actual particle size analysis of the sample material. However, the measurement results will be distorted depending on the added amount of additive.

Liquid sieving aids:

Antistatic spray, benzene, alcohol and surfactants can be used as liquid sieving aids, though benzene and alcohol are only to be used during sample preparation. They reduce the electrostatic charges, wash out fatty or oily components of the sample material, or diminish the surface tension in the wet sieving.

15 Disposal

In the case of a disposal, the respective statutory requirements must be observed. In the following, information on the disposal of electrical and electronic devices in the European Community are given.

Within the European Community the disposal of electrically operated devices is regulated by national provisions that are based on the EU Directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE).

Accordingly, all devices supplied after August 13th 2005 in the business-to-business area, to which this product is classified, may no longer be disposed of with municipal or household waste. To document this, the devices are provided with the disposal label.

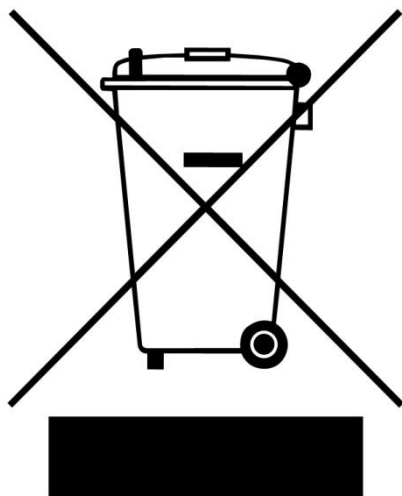


Fig. 25: Disposal label

Since the disposal regulations worldwide and also within the EU may differ from country to country, the supplier of the device should be consulted directly in case of need.

This labelling obligation is applied in Germany since March 23rd 2006. From this date on, the manufacturer must provide an adequate possibility of returning all devices delivered since August 13th 2005. For all devices delivered before August 13th 2005 the end user is responsible for the proper disposal.

16 Index

A

Accessories	50
Ambient temperature	17
Amendment status	6
Amperage	19
Amplitude	34
definition	36
dependence on the load	35
optimisation	36
optimum	36
setting	33

B

Back	31
Back view	31
Bar code	18

C

Calibration	46
Calibration service	51
Capacity	19
CE marking	18
Certificate	50
Clamping lid	30
Cleaning	47
Collecting pan with outlet	42
Complaints	16
Condensation	16
Confirmation form for the managing operator	11
Continuous operation	36
Control unit	31
connect	25
wall mounting	26
Controlling the device	33
Copyright	6

D

Degree of protection	13
Depth	14
base	14
Device designation	18
Dimensions	14
Disclaimer	6
Disconnection from the mains	22
Displays	33
Disposal	52
label	18, 52
regulations	52

E

EasySieve®	32, 44
RS232 interface	26
Electrical connection	18
Electromagnetic compatibility	14
EMC	14
Emissions	13
Equivalent continuous sound level	13

Error

E00	45
E10	45
E20	45
E26	45
Error messages	45
Evaluation	44
Explanations of the safety instructions	8
External fuse	18

F

Feed grain size	15
Feed quantity	14
First commissioning	22
Fixing nut	30
Frequency	18
Front	30
Front view	30
Functions	33
Fuse strength	19
Fuse type	19
Fuses	31
replacing	49

G

General catalogue	50
General safety instructions	9
Grain size	28
range	28
Gravitational acceleration	34

H

Hearing damage	13
Height	14
Height of lift	34
Humidity	17

I

Information note	
bS	45
ES	45
S	45
Information notes	45
Installation	16
Installation height	17
Installation site	
conditions	17
Interval	37
off	33
on	33
setting	33
Interval time	37

L

L _{eq}	13
Lifting the device	20
Load diagram	35, 43

Location requirements	14
Long-term operation	29
M	
Mains connection	31
Mains frequency	19
Mains supply	18
Mains switch	31
Maintenance	11, 46, 47, 48
Manual	6, 9, 11
Manual mode	38
Manufacturer's address	18
Materials	28
Maximum relative humidity	17
Measurement range	15
Measurement report	44
N	
Notes on the manual	6
Number of fractions maximum	15
Number of fuses	19
O	
Operating controls	33
Operating hours	39
Operating instructions	11
Operating the device	28
Organic solvents	47
P	
Packaging	16
Part number	18
Particle distribution	44
Particle size analysis	44
Particle size characteristics	44
Particle size distribution	28
Particle size range	15
Payload	15
Performing a sieving	32
Power version	19
Principle of operation	29
Process continue	34
end	34
pause	34
start	33
stop	34
Program editing	38
selecting	38
Programme saving	38
setting	33
Programme mode	38
Q	
Quick clamping lever green	30
red	31
Quick clamping unit	30

installation	25
R	
Range of application of the device	29
Rated power	14
Receptacle volume	14
Repair	10, 46, 48
Repair instructions	10
Required floor space	14
Return	16
for service and maintenance	46
Return device	52
Return form	46
RS232 interface	31
S	
Safety	7
Safety Officer	7
Sample quantity maximum	15
Scale	44
Serial number	18
Service address	10
Setting mode	32
Sieve clamping unit comfort	23
installation comfort	24
installation standard	23
standard	23
types	22
Sieve diameter	15
Sieve plate mounting	21
Sieve stack height	22
maximum height	15
maximum weight	15
Sieving aids	15, 51
Sieving noises	13
Sieving time optimum	37
Signal tone	39
Signs	6
Small accessories	50
Software	44
version	39
Sound level	13
Sound parameters	13
Spare parts	50
Spray nozzle	42
Standby mode	32
START	33
STOP	33
Switching on / off	31
Symbols	6
T	
Target group	7
Technical data	13
Temperature fluctuations	16

Temperature range	17	Vibrations	22
Temporary storage	16	Vibratory sieving	29
Test sieve	28, 50	Views of the device	30
cleaning	47	Voltage	18
diameter	22	W	
drying.....	48	Warning.....	8
maximum drying temperature	48	Information	9
selection	32	Warranty claim	16
Threaded rod	30	Warranty claims	9
Time	36	Wear.....	47, 48
optimisation	36	Wear parts	50
setting.....	33	Weight.....	14, 20
Transport	16	Wet sieving	15, 41
Transport damage	16	liquid	41
Transportation aid.....	20	perform	41
removing	19	recommended parameters	43
Type plate	18, 31	Wet sieving lid.....	42
description.....	18	Width	14
U		base.....	14
Use of the device for the intended purpose	28	Workplace related emission level	13
V		Y	
Venting ring.....	43	Year of production.....	18

VIBRATORY SIEVE SHAKER

AS 450 control | 30.026.xxxx

EU DECLARATION OF CONFORMITY

Herewith we declare, represented by the signatory, that the above mentioned device complies with the following directives and harmonized standards:

Machinery Directive 2006/42/EC

Applied standards, in particular:

DIN EN ISO 12100 Safety of machinery
DIN EN 60529 Degrees of protection provided by enclosures (IP Code)

EMC Directive 2014/30/EU

Applied standards, in particular:

DIN EN 55011 Industrial, scientific and medical equipment - Radio-frequency disturbance characteristics - Limits and methods of measurement
DIN EN 61000-3-2 Electromagnetic compatibility (EMC)
DIN EN 61000-3-3 Electromagnetic compatibility (EMC)
DIN EN 61326-1 Electrical equipment for measurement, control and laboratory use - EMC requirements

Low Voltage Directive 2014/35/EU

Applied standards, in particular:

DIN EN 61010-1 Safety requirements for electrical equipment for measurement, control and laboratory use

Authorized person for the compilation of technical documents:

Dr. Loredana Di Labio (technical documentation)

Furthermore, we declare that the relevant technical documentation for the above mentioned device has been compiled according to Annex VII Part A of the Machinery Directive, and we undertake to submit this documentation on request to the market surveillance authorities.

In case of a modification of the device not previously agreed with Retsch GmbH, as well as the use of unauthorised spare parts or accessories, this declaration will lose its validity.

Retsch GmbH



Dr. Ing. Frank Janetta, Team Leader R&D Department

Haan, 08/2017





Retsch[®]

Copyright

© Copyright by
Retsch GmbH
Retsch-Allee 1-5
42781 Haan
Germany